


# Manual

## Cryogenic Mixer Mill CryoMill



Translation

**Retsch**<sup>®</sup>

**Copyright**

© Copyright by  
Retsch GmbH  
Retsch-Allee 1-5  
42781 Haan  
Germany

## Table of Contents

<b>1</b>	<b>Notes on the Manual</b> .....	<b>7</b>
1.1	Explanations of the Safety Instructions .....	8
1.2	General Safety Instructions .....	9
1.3	Repairs.....	10
<b>2</b>	<b>Confirmation</b> .....	<b>11</b>
<b>3</b>	<b>Technical Data</b> .....	<b>12</b>
3.1	Protective Equipment.....	12
3.2	Emissions.....	12
3.3	Degree of Protection.....	13
3.4	Rated Power .....	13
3.5	Dimensions and Weight.....	13
3.6	Required Floor Space.....	13
<b>4</b>	<b>Packaging, Transport and Installation</b> .....	<b>14</b>
4.1	Packaging .....	14
4.2	Transport.....	14
4.3	Temperature fluctuations and condensed water .....	14
4.4	Conditions for the Installation Site .....	14
4.5	Installation of the Device.....	15
4.6	Removing the Transportation Lock.....	16
4.7	Type Plate Description.....	20
4.8	Electrical Connection .....	20
4.9	Connecting the Coolant Feed .....	21
4.10	Removing the liquid nitrogen feed line .....	24
4.11	Cold air outlet.....	25
4.12	Mounting the exhaust hose.....	26
<b>5</b>	<b>Operating the Device</b> .....	<b>27</b>
5.1	Use of the Device for the Intended Purpose.....	27
5.2	Operating elements and displays .....	28
5.3	Summary table of device parts .....	29
5.4	Operating Controls, Displays and Functions .....	30
5.5	Inserting the Grinding Jar .....	31
5.5.1	Inserting and removing CryoMill special grinding jar .....	32
5.5.2	Inserting the cooling jacket sealing plug .....	35
5.5.3	Insert the adapter for 4/2 x 5 ml grinding jar .....	35
5.5.4	Insert the adapter for 6 / 4 / 2 x 2-ml Eppendorf reaction vial .....	36
5.6	Explanations of the grinding cycles .....	37
5.6.1	Grinding without cooling .....	37
5.6.2	Grinding with cooling .....	38
5.7	Programme Mode .....	39
5.7.1	Grinding without program .....	39
5.7.2	Grinding with program .....	40
5.7.2.1	Selection of a program.....	40
5.7.2.2	Display of the program grinding parameters .....	40
5.7.2.3	Create/change a program.....	40
5.7.2.4	Delete program .....	41

5.8	Setting grinding cycles .....	41
5.9	Setting the precooling time .....	42
5.9.1	Precooling with definable precooling time .....	42
5.9.2	Automatic precooling .....	42
5.10	Setting the grinding time .....	43
5.11	Setting the duration of intermediate cooling .....	43
5.12	Setting the grinding frequency .....	44
5.13	Starting the grinding process .....	44
5.14	Interrupting and continuing the grinding process .....	45
5.14.1	Pausing grinding .....	45
5.14.2	Cancelling grinding .....	45
5.15	Stopping the grinding process .....	46
5.16	Operating hours display .....	46
5.17	Operating software display .....	47
5.18	Replacing the machine fuses.....	48
<b>6</b>	<b>Working instructions.....</b>	<b>49</b>
6.1	General .....	49
6.2	The grinding process .....	49
6.3	Quantities of sample material and feed sizes .....	50
6.4	Working instructions .....	50
<b>7</b>	<b>Safety functions and fault display .....</b>	<b>52</b>
7.1	Error Messages .....	52
<b>8</b>	<b>Cleaning, Wear and Maintenance .....</b>	<b>53</b>
8.1	Cleaning.....	53
8.2	Wear .....	53
8.3	Wearing parts .....	53
8.4	Maintenance .....	53
8.5	Checks.....	54
<b>9</b>	<b>Safety and occupational health instructions.....</b>	<b>55</b>
9.1	General .....	55
9.2	Safety instructions on the handling of liquid nitrogen .....	55
9.3	Oxygen deficiency .....	55
9.3.1	Dangers .....	56
9.3.2	Causes.....	56
9.3.3	Recommendations .....	56
9.3.4	General conduct in the case of an accident .....	57
9.4	Cryogenic burns.....	57
9.4.1	Dangers .....	57
9.4.2	Causes.....	57
9.4.2.1	Burns through splashes .....	57
9.4.2.2	Burns through contact.....	57
9.4.3	Recommendations .....	57
9.4.4	General rules of conduct for splashes with liquid nitrogen .....	58
9.4.4.1	To the eyes .....	58
9.4.4.2	On the skin.....	58
9.5	Danger of explosion.....	58
9.5.1	Dangers .....	58

---

9.5.2	Causes.....	58
9.5.3	Recommendations.....	58
9.5.4	General conduct in the case of an accident .....	59
9.6	Oxygen enrichment.....	59
9.6.1	Dangers .....	59
9.6.2	Causes.....	59
9.6.3	Recommendations.....	59
9.7	Surrounding area of the machine .....	60
9.7.1	Rooms.....	60
<b>10</b>	<b>Disposal.....</b>	<b>61</b>
<b>11</b>	<b>Index .....</b>	<b>62</b>



---

# 1 Notes on the Manual

**Dear user,**

please read the following manual referring to this device carefully before starting any installation, commissioning and operation.

This manual is a technical guide on how to operate the device safely and it contains all the information required for the areas specified in the table of contents. This technical documentation is a reference and instruction manual. The individual chapters are complete in themselves. Familiarity (of the respective target groups defined according to area) with the relevant chapters is a precondition for the safe and appropriate use of the device.

This manual does not contain any repair instructions. In case of any obscurities or questions with regards to this document or the device, as well as if faults arise or repairs are necessary, please contact your supplier or get in touch with Retsch GmbH directly.

Application-technical information relating to samples to be processed are not or only to a certain extend included. However, more information thereof can be found in the internet on the webpage of the respective device on the Retsch GmbH homepage (<http://www.retsch.com>).

**Revision status:**

This document revision 0004 refers to the manual "Cryogenic Mixer Mill CryoMill" in compliance with the Directive of Machinery 2006/42/EG.

**Changes**

Subject to technical changes.

**Copyright**

Disclosure or reproduction of this documentation, use and disclosure of its contents are only permitted with the express permission of Retsch GmbH.

Infringements will result in damage compensation liability.

## 1.1 Explanations of the Safety Instructions

In this Operating Manual we give you the following safety warnings

**Serious injury** may result from failing to heed these safety warnings. We give you the following warnings and corresponding content.

 **WARNING**

**Type of danger / personal injury**

Source of danger

- Possible consequences if the dangers are not observed.
- **Instructions on how the dangers are to be avoided.**

We also use the following signal word box in the text or in the instructions on action to be taken:

 **WARNING**

**Moderate or mild injury** may result from failing to heed these safety warnings. We give you the following warnings and corresponding content.

 **CAUTION**

**Type of danger / personal injury**

Source of danger

- Possible consequences if the dangers are not observed.
- **Instructions on how the dangers are to be avoided.**

We also use the following signal word box in the text or in the instructions on action to be taken:

 **CAUTION**

In the event of possible **property damage** we inform you with the word “Instructions” and the corresponding content.

*NOTICE*

**Nature of the property damage**

Source of property damage

- Possible consequences if the instructions are not observed.
- **Instructions on how the dangers are to be avoided.**

We also use the following signal word in the text or in the instructions on action to be taken:

*NOTICE*



## 1.2 General Safety Instructions


**CAUTION**
**Read the Operating Manual**

Non-observance of these operating instructions

- The non-observance of these operating instructions can result in personal injuries.
- **Read the operating manual before using the device.**
- **We use the adjacent symbol to draw attention to the necessity of knowing the contents of this operating manual.**



**Target group** : All persons concerned with the machine in any form

This machine is a modern, high performance product from Retsch GmbH and complies with the state of the art. Operational safety is given if the machine is handled for the intended purpose and attention is given to this technical documentation.

You, as the owner/managing operator of the machine, must ensure that the people entrusted with working on the machine:

- have noted and understood all the regulations regarding safety,
- are familiar before starting work with all the operating instructions and specifications for the target group relevant for them,
- have easy access always to the technical documentation for this machine,
- and that new personnel before starting work on the machine are familiarised with the safe handling of the machine and its use for its intended purpose, either by verbal instructions from a competent person and/or by means of this technical documentation.

Improper operation can result in personal injuries and material damage. You are responsible for your own safety and that of your employees.

Make sure that no unauthorised person has access to the machine.


**CAUTION**
**Changes to the machine**

- Changes to the machine may lead to personal injury.
- **Do not make any change to the machine and use spare parts and accessories that have been approved by Retsch exclusively.**

### NOTICE

**Changes to the machine**

- The conformity declared by Retsch with the European Directives will lose its validity.
- You lose all warranty claims.
- **Do not make any change to the machine and use spare parts and accessories that have been approved by Retsch exclusively.**

### 1.3 Repairs

This operating manual does not contain any repair instructions. For your own safety, repairs may only be carried out by Retsch GmbH or an authorized representative or by Retsch service engineers.

**In that case please inform:**

The Retsch representative in your country
Your supplier
Retsch GmbH directly

**Your Service Address:**

## 2 Confirmation

This operating manual contains essential instructions for operating and maintaining the device which must be strictly observed. It is essential that they be read by the operator and by the qualified staff responsible for the device before the device is commissioned. This operating manual must be available and accessible at the place of use at all times.

The user of the device herewith confirms to the managing operator (owner) that (s)he has received sufficient instructions about the operation and maintenance of the system. The user has received the operating manual, has read and taken note of its contents and consequently has all the information required for safe operation and is sufficiently familiar with the device.

As the owner/managing operator you should for your own protection have your employees confirm that they have received the instructions about the operation of the machine.

I have read and taken note of the contents of all chapters in this operating manual as well as all safety instructions and warnings.

**User**

Surname, first name (block letters)

Position in the company

Signature

**Service technician or operator**

Surname, first name (block letters)

Position in the company

Place, date and signature

### 3 Technical Data

#### 3.1 Protective Equipment

The grinding chamber of the CryoMill laboratory mill is surrounded by a stable hood. It is only possible to start the device when the hood is closed.

#### 3.2 Emissions

Noise values:

The noise values are also influenced by the properties of the sample material.

##### No-load operation:

Emission value with regard to workplace  $L_{pAeq} = 61,4 \text{ dB(A)}$

##### Example 1:

Emission value with regard to workplace  $L_{pAeq} = 75,2 \text{ dB(A)}$

##### Operating conditions:

Receptacle:	1 50ml steel grinding jar
Grinding component:	1 25mm steel ball each
Feed material:	Quartz gravel particles approx. 4.0 – 6.0mm
Feed quantity:	8ml
Frequency:	30 Hz

##### Example 2:

Emission value with regard to workplace  $L_{pAeq} = 66,3 \text{ dB(A)}$

##### Operating conditions:

Receptacle:	4 5ml steel grinding jar
Grinding component:	2 8mm tungsten carbide balls each
Feed material:	Quartz gravel particles approx. 1.0 – 1.5mm
Feed quantity:	1.5ml
Frequency:	30 Hz



**Hearing loss**

A high sound level may occur depending on the type of material, the number of balls used, the set grinding frequency and the grinding time.

- Excess noise in terms of intensity and duration can lead to impairments or permanent damage to hearing.
- **Ensure you take suitable soundproofing measures or wear hearing protection.**



**3.3 Degree of Protection**

IP30

**3.4 Rated Power**

260 watts

**3.5 Dimensions and Weight**

Height : 373 to approx. 630 mm with open hood

Width: 395 mm

Depth: 577 mm/ 710mm with exhaust hose

Weight: approx. 46 kg without grinding jar

**3.6 Required Floor Space**

Installation surface: 620 mm x 720 mm

## 4 Packaging, Transport and Installation

### 4.1 Packaging

The packaging is designed for the transport route. It corresponds to generally acceptable packaging guidelines.

---

#### *NOTICE*

##### **Storage of packaging**

- In the event of a complaint or return, your warranty claims may be endangered if the packaging is inadequate or the machine has not been secured correctly.
  - **Please keep the packaging for the duration of the warranty period.**
- 

### 4.2 Transport

---

#### *NOTICE*

##### **Transport**

- Mechanical or electronic components may be damaged.
  - **The machine may not be knocked, shaken or thrown during transport.**
- 

### 4.3 Temperature fluctuations and condensed water

---

#### *NOTICE*

##### **Temperature fluctuations**

The machine may be subject to strong temperature fluctuations during transport (e.g. aircraft transport)

- The resultant condensed water may damage electronic components.
  - **Protect the machine from condensed water.**
- 

### 4.4 Conditions for the Installation Site

Ambient temperature: 5°C to 40°C

---

#### *NOTICE*

##### **Ambient temperature**

- Electronic and mechanical components may be damaged and the performance data alter to an unknown extent.
  - **Do not exceed or fall below the permitted temperature range of the machine (5°C to 40°C / ambient temperature).**
-

Atmospheric humidity:

Maximum relative humidity 80% at temperatures up to 31°C, decreasing linearly up to 50% relative humidity at 40°C

---

*NOTICE*

**Atmospheric humidity**

- Electronic and mechanical components may be damaged and the performance data alter to an unknown extent.
  - **Do not exceed the admissible range for atmospheric humidity.**
- 

## 4.5 Installation of the Device

Installation height: maximum 2000 m above sea level

---

*NOTICE*

**Property damage**

A very cold vapour plume escapes from the nitrogen gas outlet.

- The vapour plume can greatly cool objects. This is why the humidity in the air condenses.
  - **The area above the nitrogen gas outlet (B) up to the ceiling must be kept free.**
- 

---

*NOTICE*

**Installation**

- Depending on the operating status of the mill, there may be slight vibrations.
  - **Place the mill on an even, flat and balanced supporting surface only. The supporting surface must be stable and must not vibrate.**
- 

---

*NOTICE*

**Installation of the machine**

- It must be possible to disconnect the machine from the mains at any time.
  - **Install the machine such that the connection for the mains cable is easily accessible.**
-

## 4.6 Removing the Transportation Lock

*NOTICE*

N1.0018

**Transportation lock**

Transport without transportation lock, or operation with transportation lock

- Mechanical components may be damaged.
- **Only transport the device with mounted transportation lock.**
- **Do not operate the device with built-in transportation lock.**

- Place the device on a stable table.
- Pull the device carefully beyond the edge of the table in order to be able to reach the screw (TA).

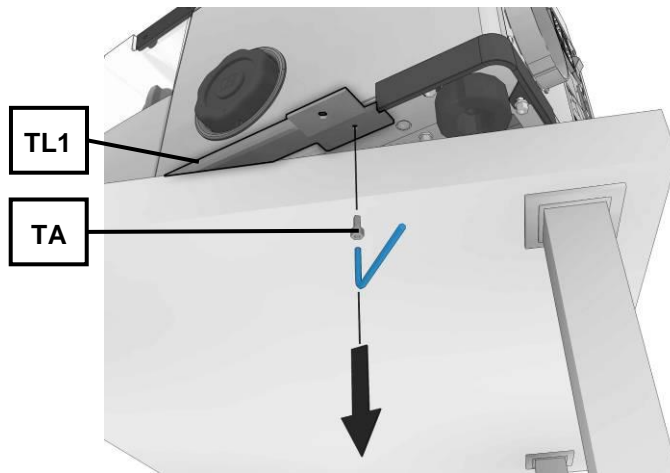


Fig. 1: Removing the transport bracket

- Remove the screw (TA) using an Allen key.

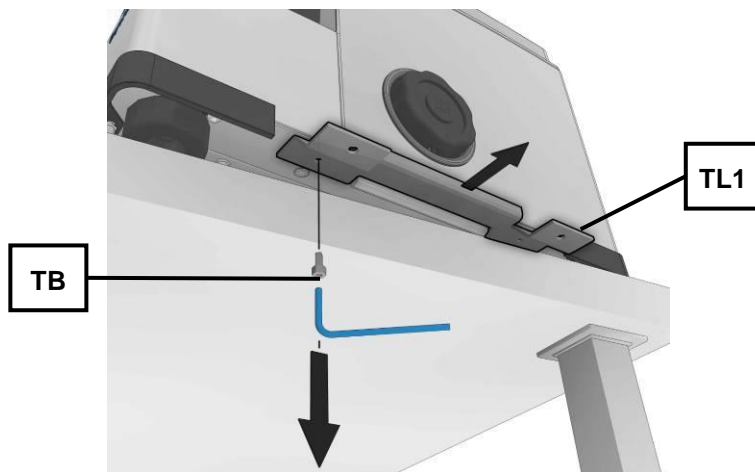


Fig. 2

- Pull the device into the position shown in order to be able to remove the screw (TB).



- Remove the second screw (**TB**) and remove the transport bracket (**TL1**).

*NOTICE*

Keep the transport brackets, screws and transport lock for future transport!

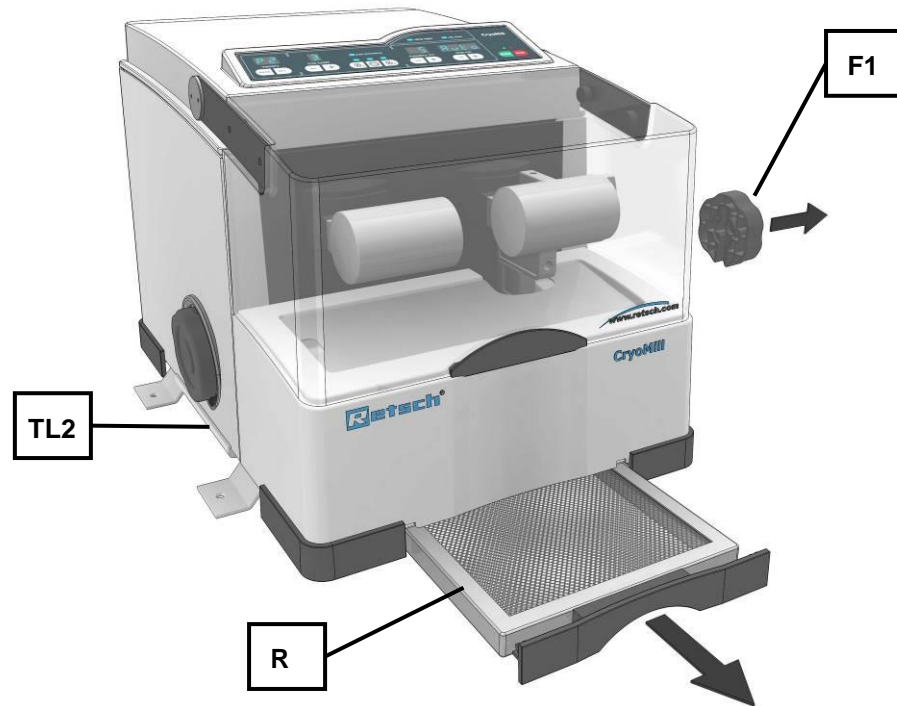


Fig. 3: Pulling out the condensation collecting filter

- Pull the condensation collecting filter (**R**) completely out of the device.
- Remove the opening aids (**F1**) on the side on which you have removed the transport bracket (**TL1**).
- Place a soft, clean mat (**U**) next to the device.
- Tilt the device onto the side on which the transport bracket (**TL1**) has been removed.

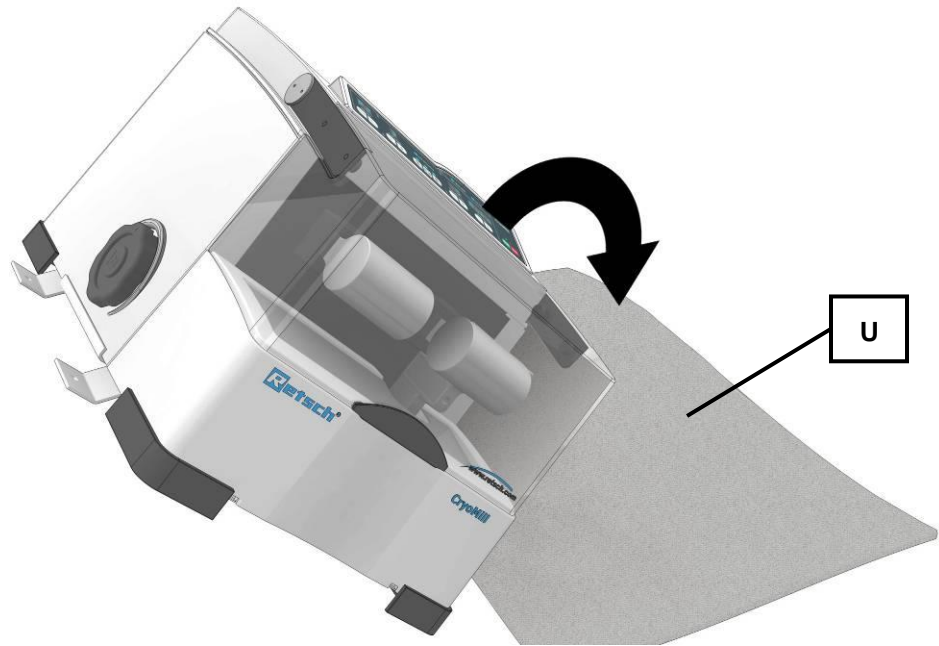


Fig. 4: Tilting the device onto its side

The transport lock (**TS**), which is labelled with an arrow, is situated on the underneath of the device.

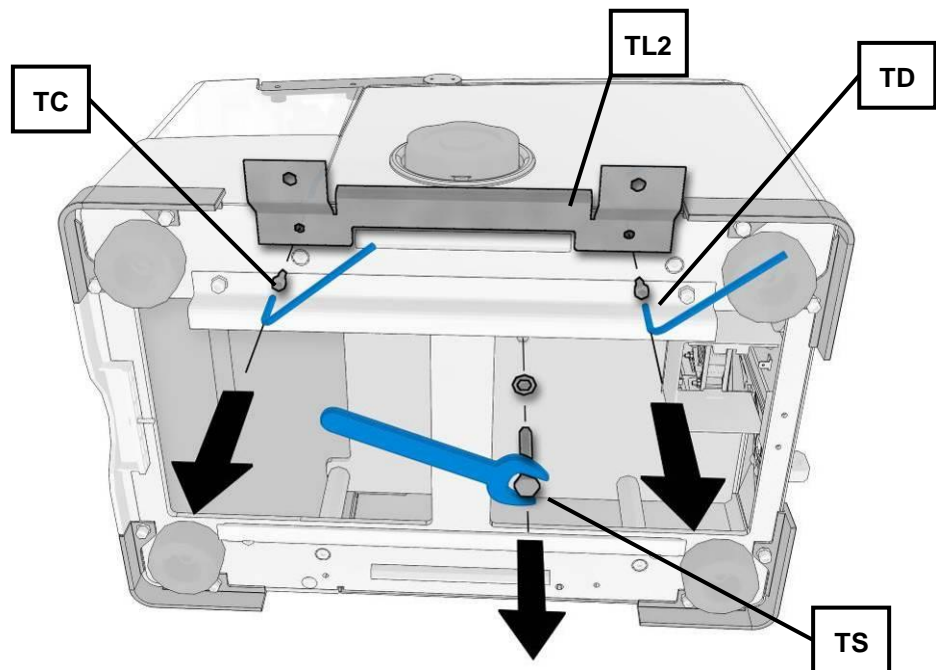


Fig. 5: Removing the transport lock

- Loosen and remove the screw (**TS**).
- Unscrew the two screws (**TC**) and (**TD**).
- Remove the transport bracket (**TL2**).

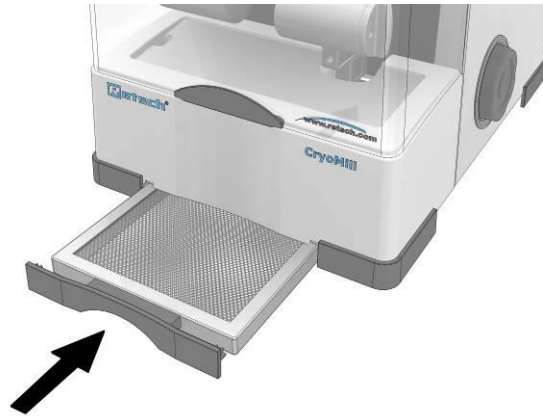


Fig. 6

- Re-insert the condensation collecting filter (R).

*NOTICE*

Keep the transport brackets, screws and transport lock for future transport!

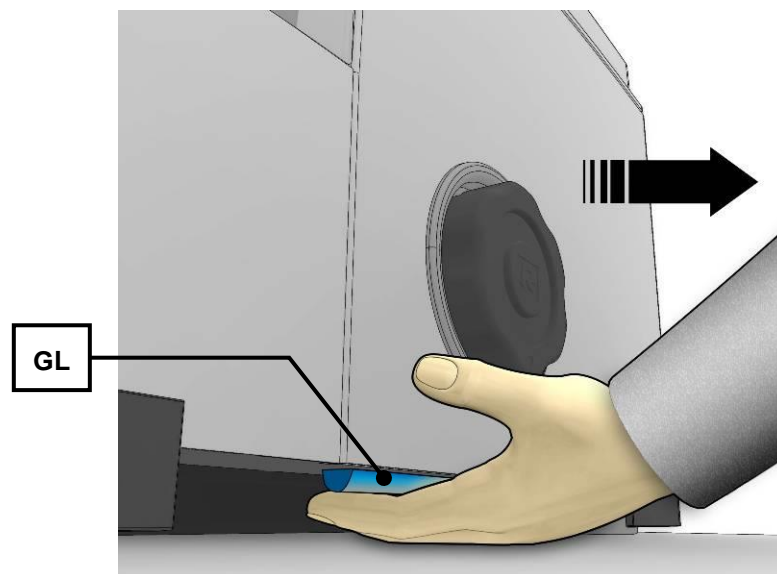


Fig. 7: Grip

A grip (GL) is situated underneath both sides of the device

- Use the grip (GL) to move the device.

### 4.7 Type Plate Description

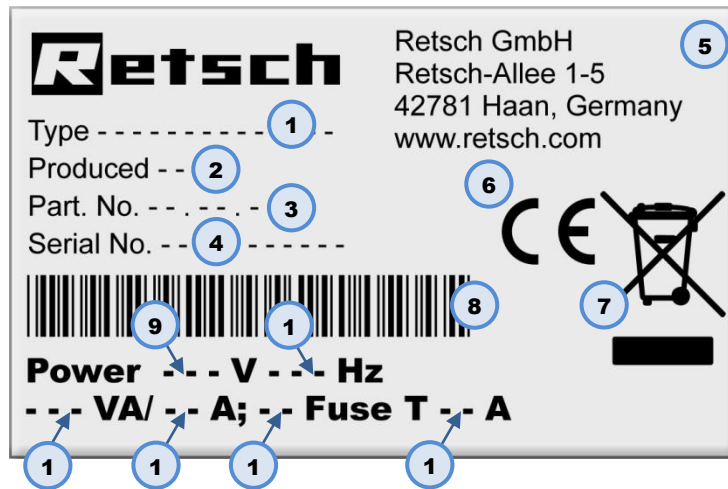


Fig. 8: Type plate lettering

- 1 Device designation
- 2 Year of production
- 3 Part number
- 4 Serial number
- 5 Manufacturer's address
- 6 CE marking
- 7 Disposal label
- 8 Bar code
- 9 Power version
- 10 Mains frequency
- 11 Capacity
- 12 Amperage
- 13 Number of fuses
- 14 Fuse type and fuse strength

In the case of questions please provide the device designation (1) or the part number (3) and the serial number (4) of the device.

### 4.8 Electrical Connection

**⚠ WARNING**

When connecting the power cable to the mains supply, use an external fuse that complies with the regulations applicable to the place of installation .

- Please check the type plate for details on the necessary voltage and frequency for the device.
- Make sure the levels agree with the existing mains power supply.
- Use the supplied connection cable to connect the device to the mains power supply.

The external fuse must be at least T6.3A (230V) T8A (100/120V).

 **WARNING**

W0002

**Danger to life through electric shock**

- An electric shock can cause burns and cardiac arrhythmia, as well as respiratory failure and cardiac arrest.
- **Never use a damaged power cable to connect the device to the mains.**
- **Check the power cable and plug for any damage before use**


*NOTICE***Electrical connection**

- Mechanical or electronic components may be damaged.
- **Please observe the information on the type plate.**

**4.9 Connecting the Coolant Feed**
 **CAUTION**
**Blocked pipe**

## Ice formation

- Water vapour from the environment may condense on parts of the device when the climate in the room or area has been altered due to transportation or changes in the air conditioning. This applies in particular following transportation of the device by plane.
- Ice or foreign bodies may block the pipes and then lead to malfunctions.
- **Keep the cooling system and inlet hose dry and free from foreign bodies.**
- **Do not use the grinding jar if it is damp or covered with frost.**
- **Place the cooling jacket seal (KA) if the device is not going to be used for a longer period of time.**
- **Keep climatic conditions in the room where the device is operated as constant as possible.**
- **Allow a 24-hour acclimatisation period following any changes to the climatic conditions for the device.**
- **The device should not be operated unsupervised.**

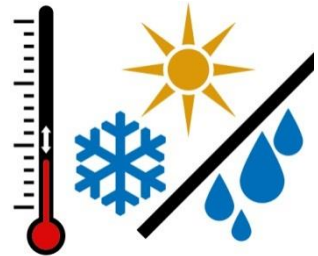

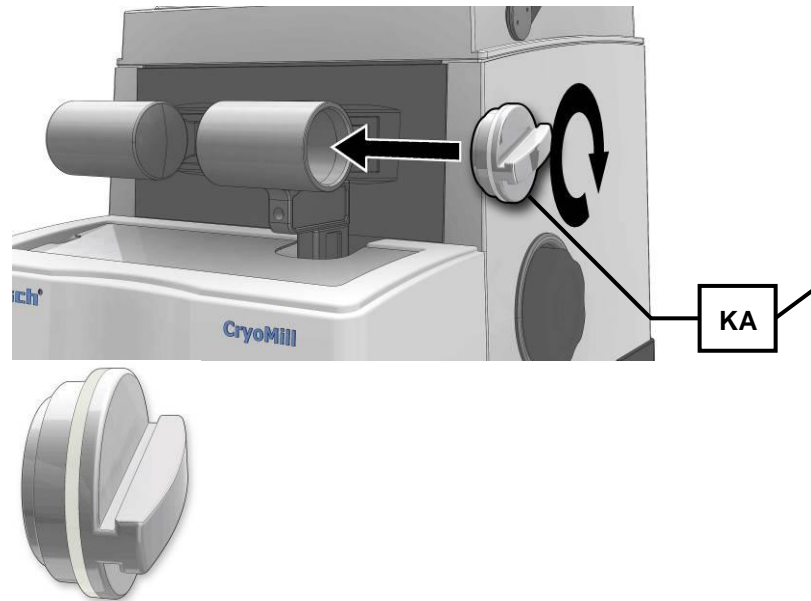


Fig. 9: Cooling jacket seal

**NOTICE**

**Damage to the machine**

Incorrect cooling agent

- The use of cooling agents other than liquid nitrogen will cause the cooling system to lose its leak tightness.
- **Only liquid nitrogen (LN2) may be used as cooling agent with this device.**

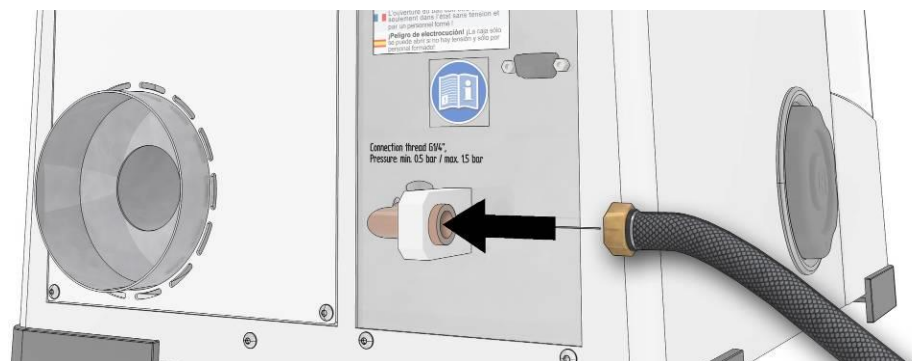


Fig. 10: Connecting the cooling agent supply

**NOTICE** After refilling the liquid nitrogen container, wait at least 12 hours (until the hose is completely dry) before connecting the container to the CryoMill. Condensed air humidity and ice crystals inside the hose might block the valves.

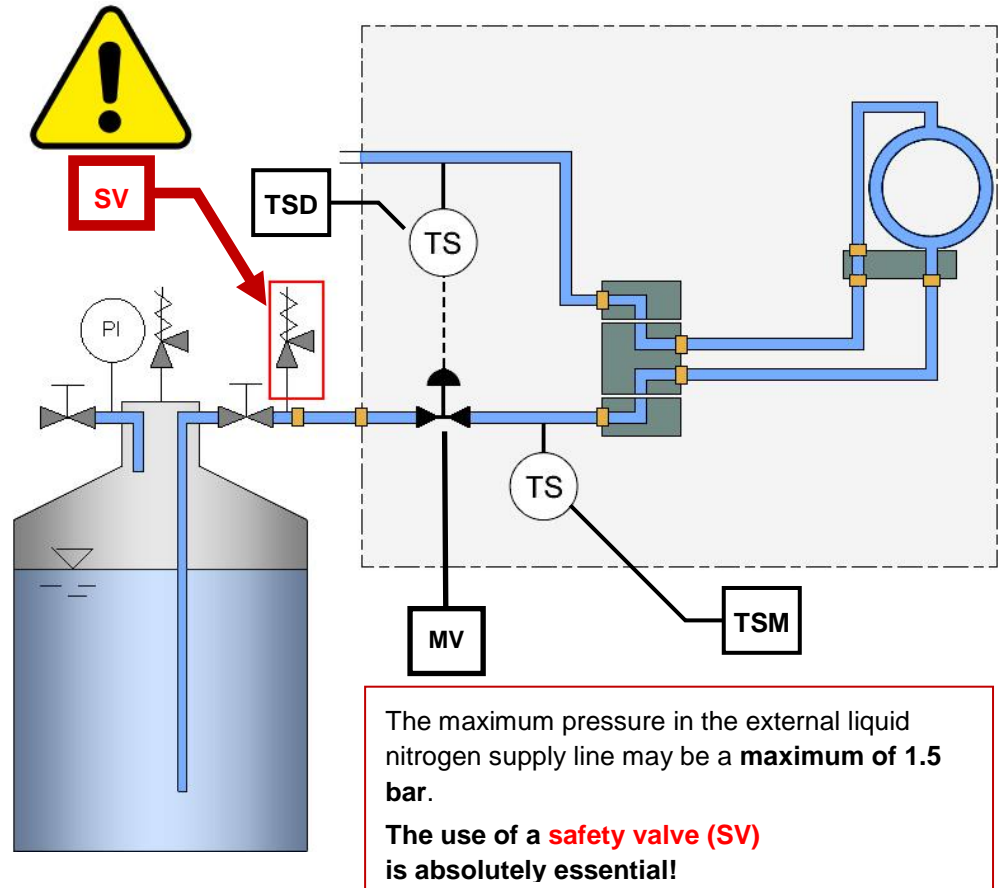


Fig. 11: Diagram of the cooling system

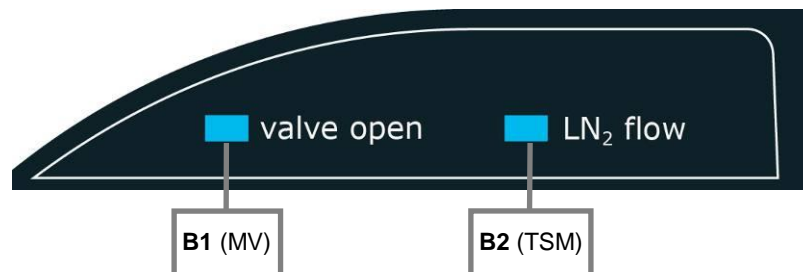


Fig. 12: LED displays

The temperature sensor (**TSD**) in the LN<sub>2</sub> outlet controls the solenoid valve (**MV**). The LED (**B1**) lights up when the solenoid valve (**MV**) is open.

The temperature sensor (**TSM**), which is downstream of the solenoid valve, registers whether liquid nitrogen is flowing into the CryoMill. The LED (**B2**) lights up when liquid nitrogen flows into the cooling system.

**CAUTION**

The device must not be moved in a cryogenic state and no additional forces may act on the cooling agent connection (**N**). The cooling agent inlet hose must not be moved while the device is in a cryogenic state.

 **CAUTION**

**Risk of injury and danger of frostbite**

Bursting of the inlet hose.

- Depending on the operating state of the device, liquid nitrogen may be enclosed in the supply line. When heated this may lead to the inlet hose bursting.
- **The maximum pressure in the external liquid nitrogen supply line may be a maximum of 1.5 bar. The use of a safety valve is absolutely essential!**

The connection for the liquid nitrogen (LN2) is a G 1/4 inch threaded connector .

The operating pressure in the liquid nitrogen supply line should be at least 0.5 bar. The lower the pressure of the liquid nitrogen supply line, the longer the precooling takes. At 0.5 bar supply line pressure and with the system at room temperature, precooling requires approx. 10 minutes. Depending on the application, you can execute up to 4 grinding operations with 10 litres of liquid nitrogen.

Throughout the entire duration of cooling, the air humidity in the chamber condenses on all cooled parts as a layer of ice .

After cooling this condensation drips into the collecting filter and evaporates.

 **CAUTION**

1.V0070

**Risk of injury and danger of frostbite**

Uncontrolled escape of liquid nitrogen

- There is a risk of injury in the event of an uncontrolled escape of liquid nitrogen.
- **Stop the liquid nitrogen supply line immediately!**

*NOTICE*

Depending on the ambient temperature, air humidity and duration of use, small quantities of water may drip from the collecting filter.

**4.10 Removing the liquid nitrogen feed line**

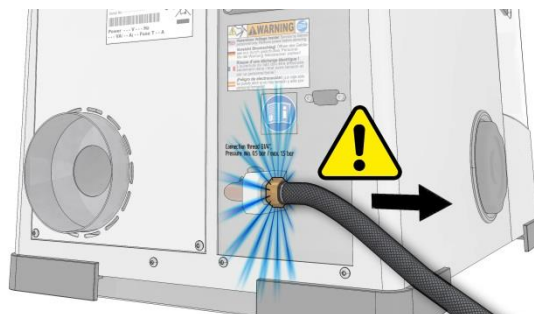


Abb. 13: Removing the liquid nitrogen feed line



**CAUTION**

Before you remove the liquid nitrogen feed line, the pressure must be relieved for safety reasons.

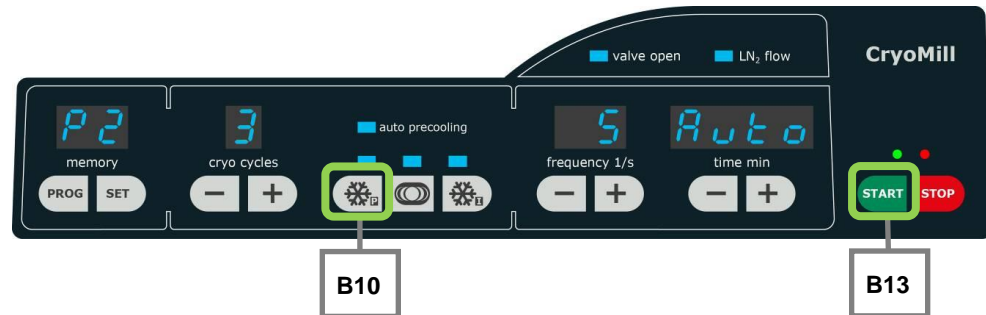


Abb. 14: pressure relief

- Press and hold the pre-cooling button (**B10**) directly followed by the START button (**B13**). Hold down both buttons for two seconds.
- As long as you keep the keys depressed, the solenoid valve will remain open.

#### 4.11 Cold air outlet

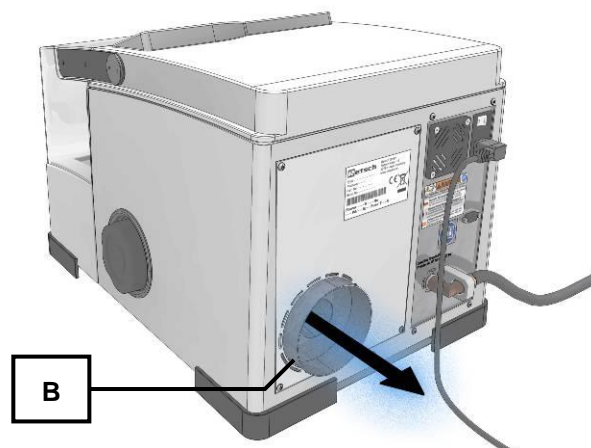


Fig. 15: Outlet for gaseous nitrogen (**B**)

**CAUTION**

**Risk of injury to eyes and skin**

Outlet for gaseous nitrogen

- Liquid nitrogen has a temperature of  $-196\text{ °C}$  and can cause burn-like injuries or frostbite in the case of contact with the skin or eyes.
- **Do not keep any parts of the body in front of the outlet for gaseous nitrogen (**B**).**

#### 4.12 Mounting the exhaust hose

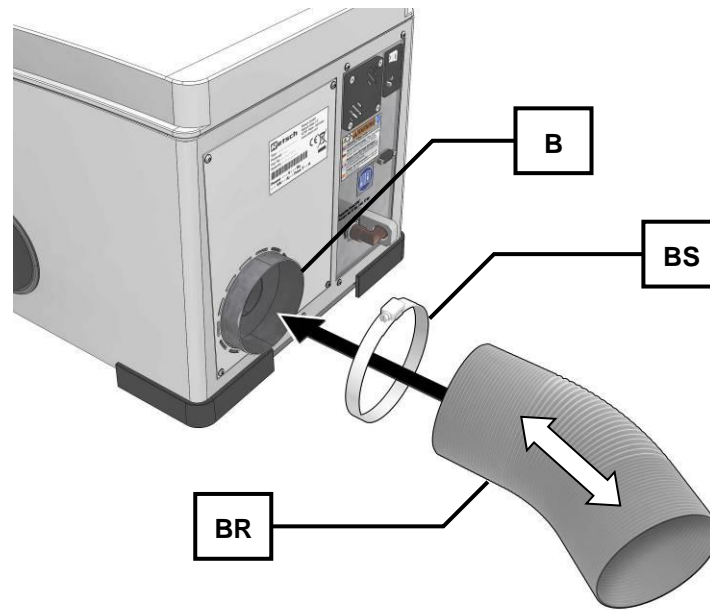


Fig. 16: Exhaust hose

- Fasten the exhaust hose (**BR**) with the clamp (**BS**) to the outlet (**B**).
- The exhaust hose can be extended up to 1 m

---

## 5 Operating the Device

### 5.1 Use of the Device for the Intended Purpose

**Target group:** operating company, operator

**Name of machine model:** CryoMill

This machine is a laboratory device and is suitable for the particularly gentle grinding and homogenisation of heat-sensitive, soft, fibrous, hard and brittle materials in a dry and wet state.

The CryoMill is designed exclusively for special grinding jars with screw-top lids (no standard grinding jars) with steel jacket. It is used for the fast pulverization of up to six samples simultaneously (depending on the type of grinding jar).

Available grinding jar volumes :

- 50ml;
- 35ml;
- 25ml;
- 10ml;
- 4 x 5ml including adapter; and
- 6 x 2ml Eppendorf vials.

The closed grinding system guarantees complete recovery of the samples . Due to the extremely short grinding time and the high final fineness of the grinding sample, the CryoMill is also ideally suitable for sample preparation for all spectral analyses.

Final fineness levels of up to 5 µm can be achieved depending on the grinding time and the specific properties of the sample material.

The optimal grinding jar filling is generally 1/3 of the grinding jar volume. Voluminous materials such as wool, leaves, grasses and similar are exceptions to this. A fill level of 70 - 80% is necessary here.

---

#### *NOTICE*

##### **Area of use of the machine**

- This machine is a laboratory machine designed for 8-hour single-shift operation.
  - **This machine may not be used as a production machine nor is it intended for continuous operation.**
-

5.2 Operating elements and displays

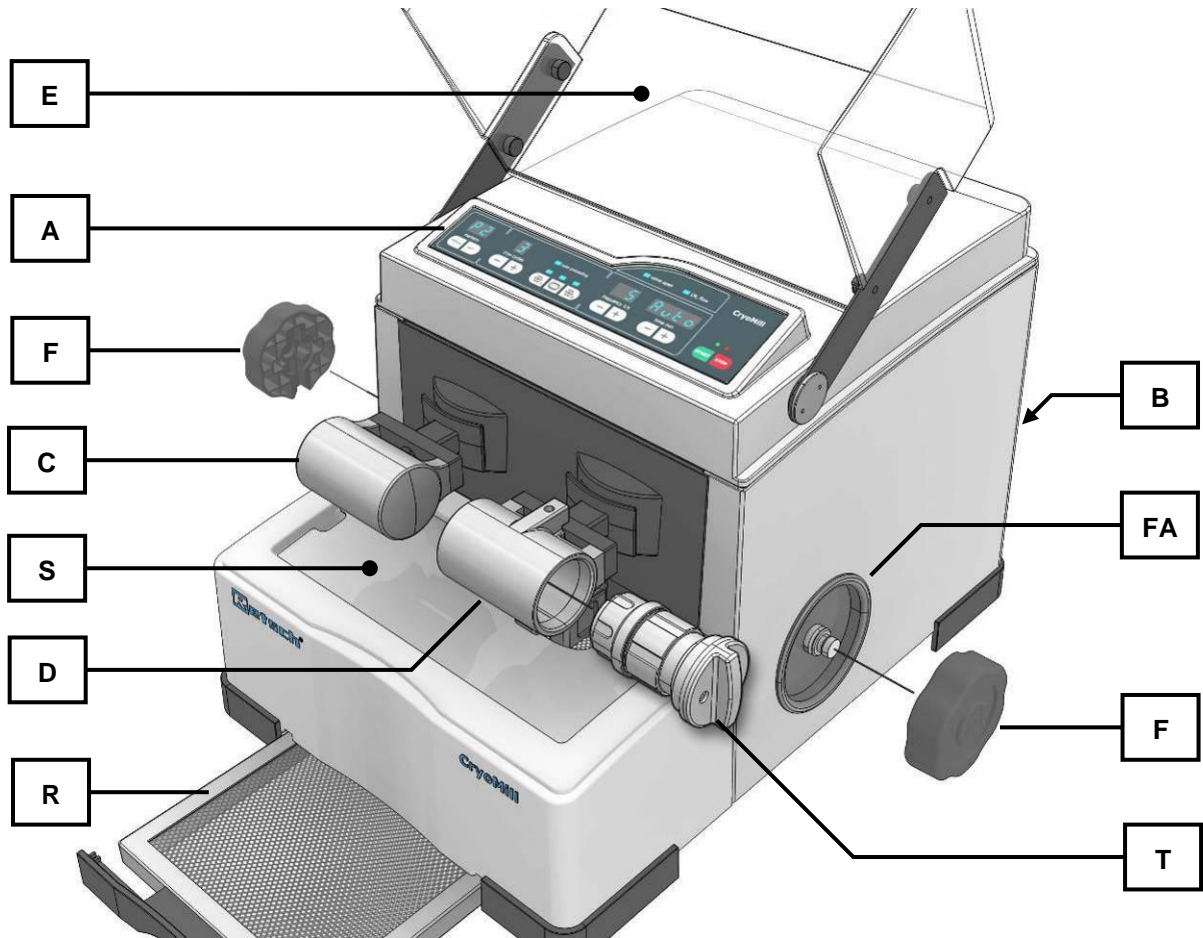


Fig. 17: Front view

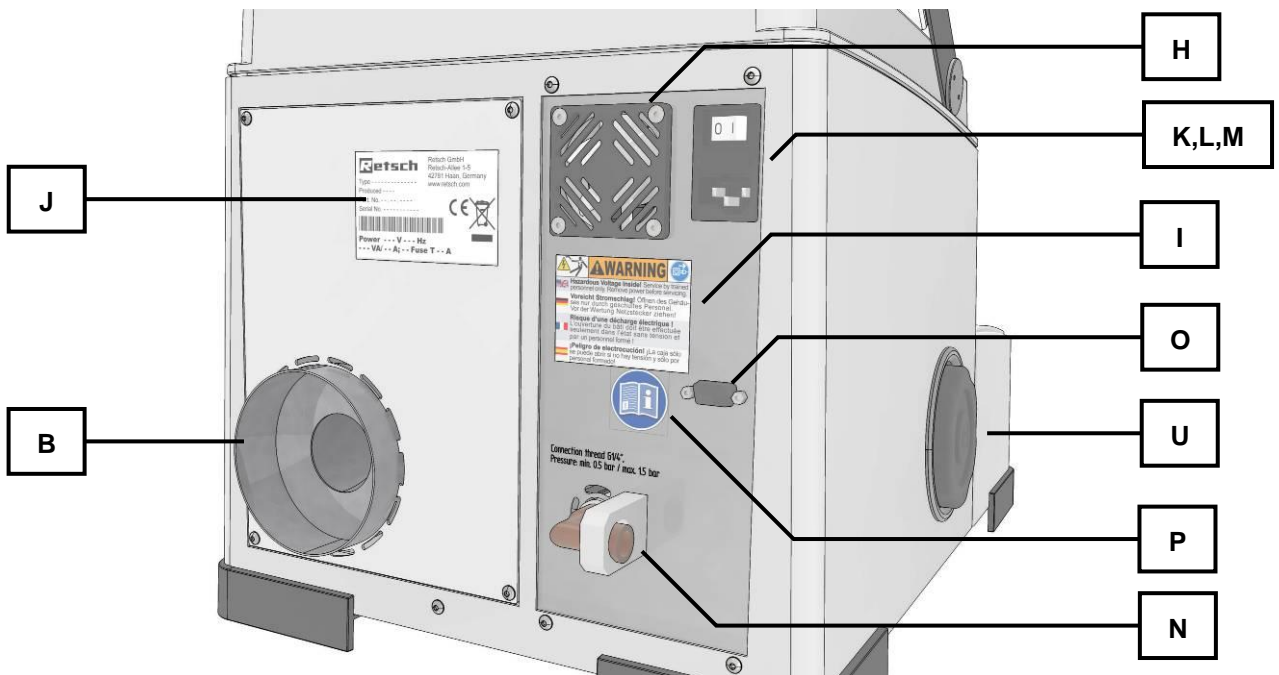


Fig. 18: Rear view

### 5.3 Summary table of device parts

Element	Description	Function
<b>A</b>	Display and control unit: See below for explanations	Time preselection, frequency preselection and starting/stopping the machine
<b>B</b>	Nitrogen gas outlet	Outlet for gaseous nitrogen
<b>C</b>	Counterweight left	Counterweight for the cooling jacket with grinding jar
<b>D</b>	Cooling jacket	Forms the cooling system together with the grinding jar ( <b>T</b> )
<b>E</b>	Hood	Closes the grinding chamber
<b>F</b>	Opening aid	Can be mounted on the grinding jar cover for simple removal of the grinding jar
<b>FA</b>	Support for opening aid	Storage of the opening aid
<b>H</b>	Fan	Ventilates the motor and the interior of the mill
<b>I</b>	Caution remove mains plug sign	Safety instruction
<b>J</b>	Type plate	Information about the device and connected loads
<b>K</b>	ON / OFF switch	Switches the device on and off
<b>L</b>	Fuse drawer	Accommodates two glass fuses
<b>M</b>	Connector	Connector for the mains lead of the device
<b>N</b>	Cooling agent connection	Connection for liquid nitrogen
<b>O</b>	RS232 interface	Enables updating of the operating software
<b>P</b>	Observe operating instructions sign	Safety instruction
<b>R</b>	Condensation collecting filter	Collects the condensation from the cooling system and allows it to evaporate
<b>S</b>	Drip tray	Housing protection
<b>T</b>	Grinding jar	Forms the cooling system together with the cooling jacket ( <b>D</b> )

5.4 Operating Controls, Displays and Functions

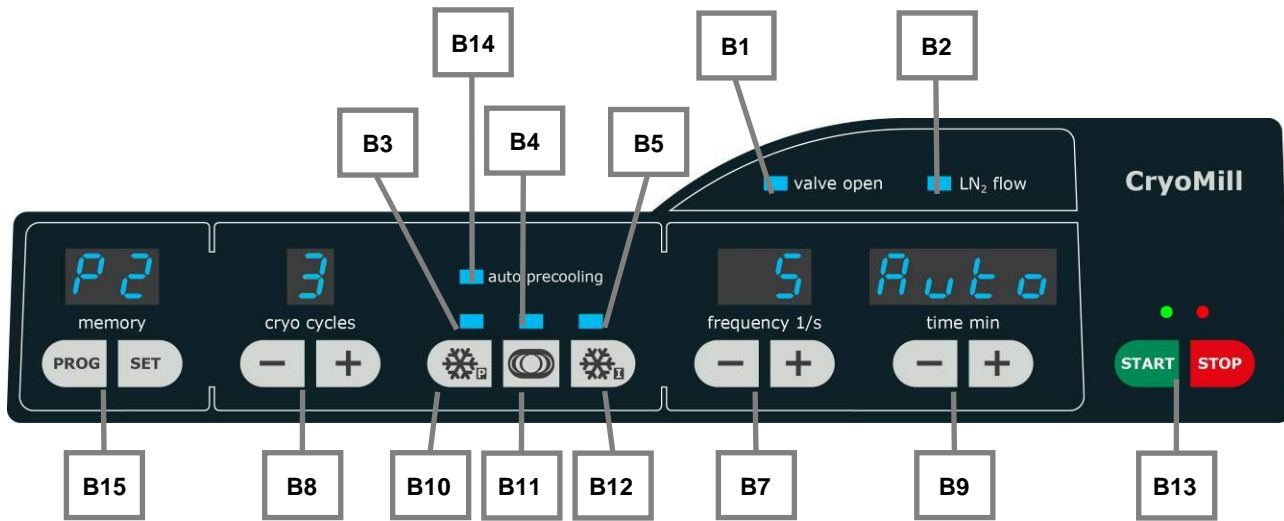


Fig. 19: Control unit

Element	Description	Function
B1	LED on – valve open LED off – valve closed	Indicates whether the internal solenoid valve on the liquid nitrogen supply line is open or closed.
B2	LED on – liquid nitrogen is flowing LED off – no liquid nitrogen	Indicates that liquid nitrogen is flowing into the device.
B3	LED precooling is running (Vkz)	Precooling is taking place / lights up when values are set
B4	LED grinding is running (Mz)	The sample is ground / lights up when values are set
B5	LED intermediate cooling is running (Zkz)	Intermediate cooling is taking place / lights up when values are set
B7	Frequency minus button	In setting mode and in operating mode: reduces the set frequency by 1 Hz (1/s). Continuous pressing switches on fast scrolling of numbers.
	Frequency plus button	In setting mode and in operating mode: increases the set frequency by 1 Hz (1/s). Continuous pressing switches on fast scrolling of numbers.
B8	Cooling cycle minus button	Reduces the number of total cycles by 1
	Cooling cycle plus button	Increases the number of total cycles by 1
B9	Time minus button	In setting mode and in operating mode: Reducing the set <ul style="list-style-type: none"> <li>• Precooling</li> <li>• Grinding time</li> <li>• Intermediate cooling time</li> </ul>

	Time plus button	In setting mode and in operating mode: Increasing the set <ul style="list-style-type: none"> <li>• Precooling</li> <li>• Grinding time</li> <li>• Intermediate cooling time</li> </ul>
<b>B10</b>	Setting precooling	Preselection button to set the precooling / switching automatic precooling on or off)
<b>B11</b>	Setting grinding	Preselection button to set the grinding time
<b>B12</b>	Setting intermediate cooling	Preselection button to set the intermediate cooling time
<b>B13</b>	Start button Green LED	Starts grinding operation Indicates grinding operation
	Stop button Red LED and ON	Interrupts or ends grinding operation, places the machine in stand-by mode
<b>B14</b>	LED - auto precooling	Indicates whether the automatic precooling is switched on
<b>B15</b>	PROG button	Button to select a program
	SET button	Button starts / saves a program editing

## 5.5 Inserting the Grinding Jar



### Danger of injury and frostbite

Escape of liquid nitrogen

- Liquid nitrogen has a temperature of  $-196\text{ }^{\circ}\text{C}$  and may cause injuries similar to burns on skin or eye contact or frostbite.
- **Never switch the machine on without tightly closed cooling casing.**



### Crushed or bruised fingers

Falling grinding chamber protective hood

- The protective hood of the grinding chamber can cause crushed or bruised fingers if it falls down.
- **Hold the flap tight when closing.**

### 5.5.1 Inserting and removing CryoMill special grinding jar

 **CAUTION**

**Clogged piping**

Formation of ice

- Ice or foreign matter can clog the pipe and subsequently cause malfunctioning.

**Keep the cooling system and the supply pipe dry and free of foreign matter.**

- **Do not insert the grinding jar if it is damp or covered with frost.**
- **Use the cooling jacket sealing plug (KA) if you are not going to use the device for a long time.**
- **Do not let the device run without supervision.**

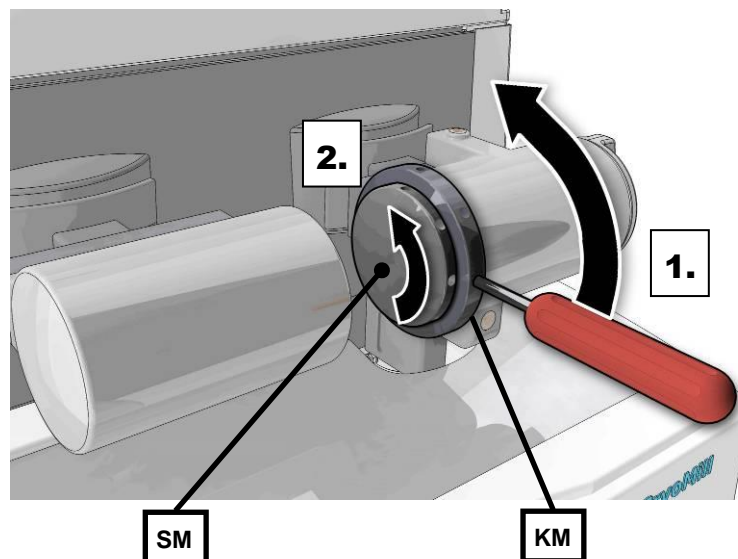


Fig. 20: Dismantling the grinding jar retaining device

- Unscrew the lock nut (KM).
- Remove the screw (SM).



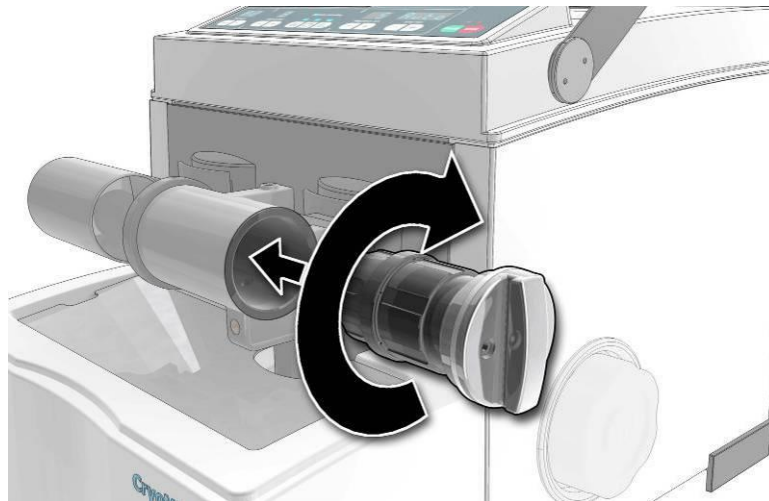


Fig. 21: Inserting the grinding jar

- Insert the special grinding jar filled with the grinding material and grinding balls into the cooling jacket .
- Twist the grinding jar in completely.

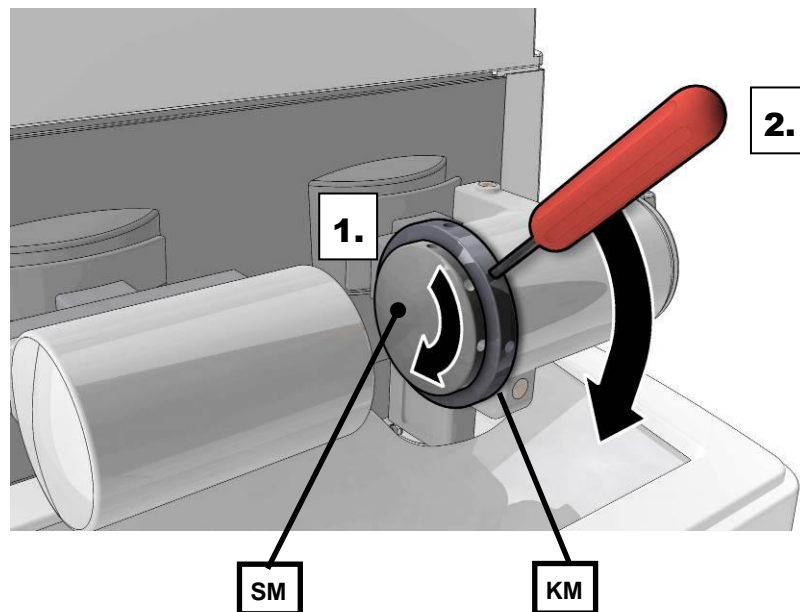


Fig. 22: Securing the grinding jar

- First tighten the screw (SM) by hand (hand-tight).
- Then tighten the lock nut against the cooling jacket. Use the locking pin provided to do this. Ensure that the grinding jar is tightly screwed to the cooling jacket, otherwise liquid nitrogen can escape.

The lock screw (SM) must be screwed tight and secured with the lock nut so that no sample material escapes from the grinding jar.

 **CAUTION**

**Risk of injury to eyes and skin**

The cooling jacket and grinding jar reach very low temperatures during grinding.

- Risk of injury to eyes and skin from severe frostbite
- **As a matter of principle, always wear goggles and protective gloves when opening the cooling jacket and the grinding jar.**



Use the opening aid (F) that is kept on the side of the device to remove the grinding jar .

- Unscrew the lock nut (KM)
- Remove the screw (SM)
- Mount the opening aid (F) on the grinding jar and open it.

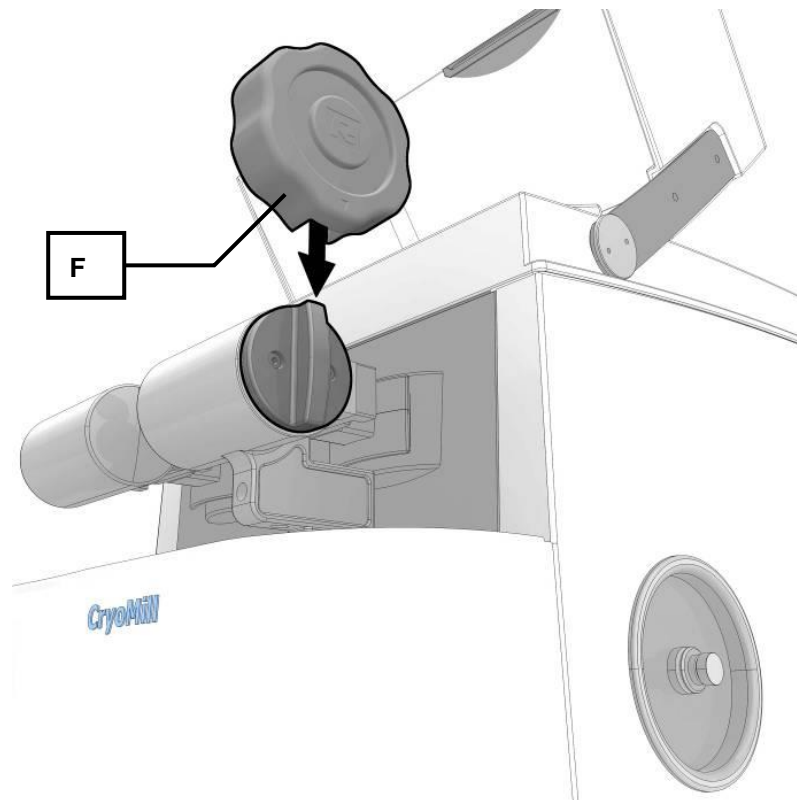


Fig. 23: Using the opening aid

**NOTE**

H0017

**Damage to the PTFE grinding jar**

- The threading on the grinding jar can become damaged.
- **When it is being inserted into the device, the PTFE grinding jar must be at room temperature and must not be pre-cooled.**

### 5.5.2 Inserting the cooling jacket sealing plug

**NOTICE**

**Blocked pipe**

Ice formation

- Ice or impurities can block pipes and thereby cause malfunctions.
- **Insert the cooling jacket seal (KA) as soon as you are no longer using the device.**

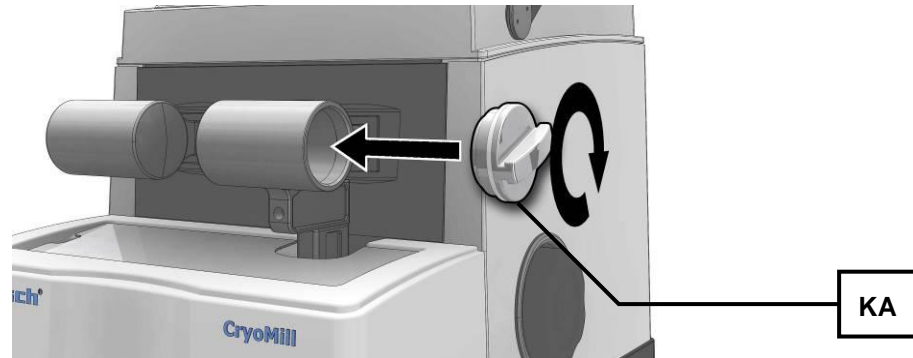


Fig. 24: Inserting the cooling jacket seal

### 5.5.3 Insert the adapter for 4/2 x 5 ml grinding jar

You can operate the CryoMill with 4 or 2 (5 ml) grinding jars.

- Insert the 5 ml grinding jars into the adapter as shown in the following diagram.

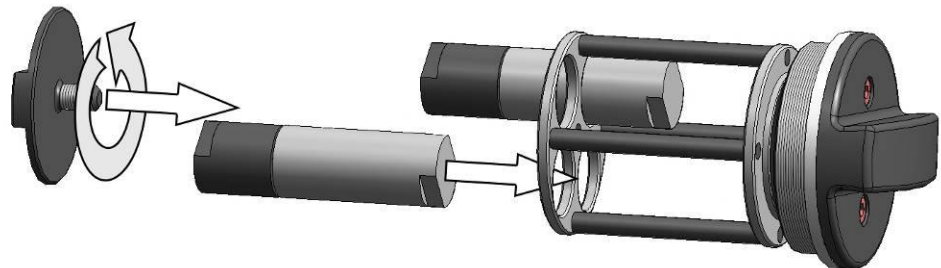


Fig. 25: Inserting a 5ml grinding jar

Note the positioning of the grinding jars. If you are using two grinding jars, these must be inserted in opposite openings.

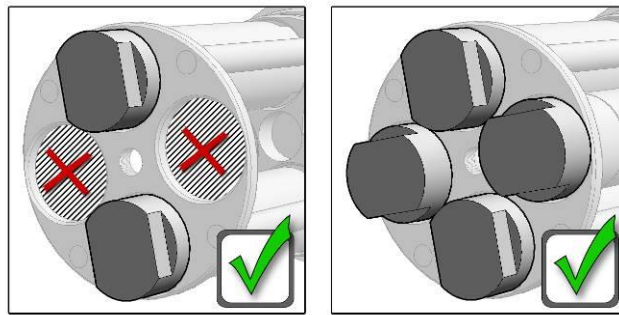


Fig. 26: 5ml grinding jar position

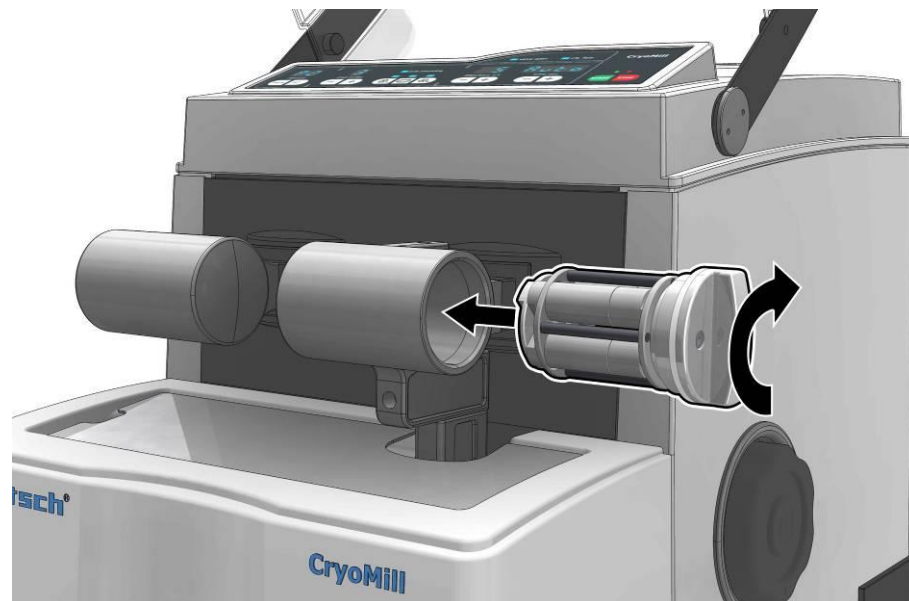


Fig. 27: Inserting the grinding jar adapter

#### 5.5.4 Insert the adapter for 6 / 4 / 2 x 2-ml Eppendorf reaction vial

##### NOTICE

##### Breaking of Eppendorf reaction vials

##### Changing material properties

- The Eppendorf reaction vials that are available as accessories change their material properties during the cryogen grinding.
- **Do not reuse reaction vials after the cryogen grinding.**
- **Do not insert reaction vials in centrifuges after the cryogen grinding**
- **Grinding may take place with a maximum 25Hz.**

You can operate the CryoMill with 2, 4 or 6 Eppendorf reaction vials (2ml).

**NOTICE**

Grinding with Eppendorf reaction vials may only take place with a maximum 25Hz. (**frequency 1/s**)

- Insert the 2ml Eppendorf reaction vials into the adapter for reaction vials as shown in the diagram below.

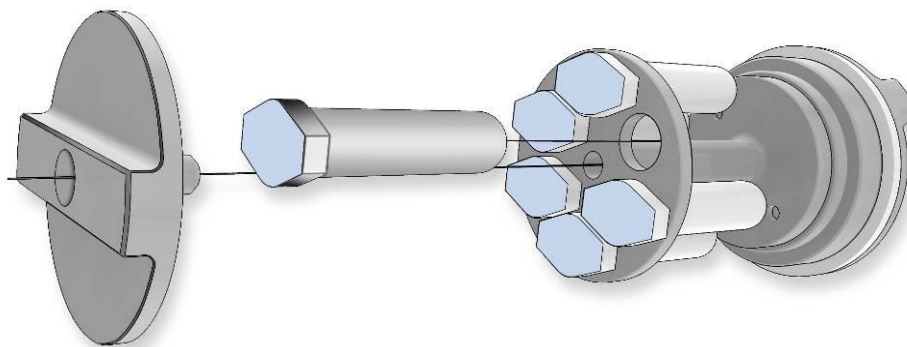


Fig. 28: Inserting 2ml Eppendorf reaction vials

Note the positioning of the Eppendorf reaction vials. If you are using two or four reaction vials, these must be inserted in opposite openings.

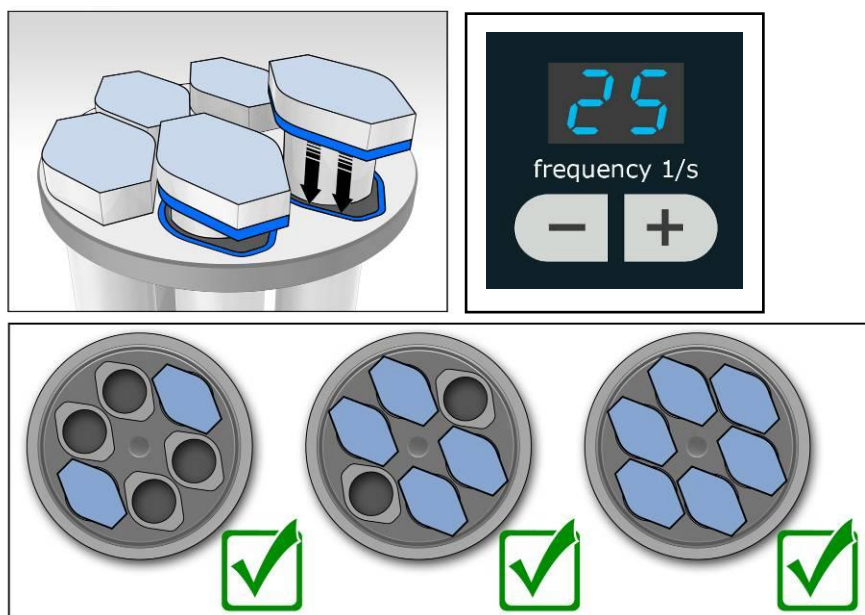


Fig. 29: Position of 2 ml reaction vials

## 5.6 Explanations of the grinding cycles

You can deploy the CryoMill for grinding with cooling or grinding without cooling.

### 5.6.1 Grinding without cooling

- Switch the CryoMill on at the main switch.
- Set the cryo cycles to - .



- Set the desired grinding time (B9)+(B11).
- Press the START button.

The START LED lights up. The remaining grinding time and the set frequency are displayed.

### 5.6.2 Grinding with cooling

The following program points can be set when grinding with cooling.

- precooling (Vkz)
- grinding time (Mz)
- intermediate cooling time (Zkz)
- cryo cycles (Number of cooling/grinding cycles)
- frequency 1/s (Grinding frequency)

A single grinding cycle consists of the precooling and the set grinding time. The intermediate cooling time does not apply in the case of a single grinding cycle.

$$Vkz + Mz = \text{total grinding time}$$

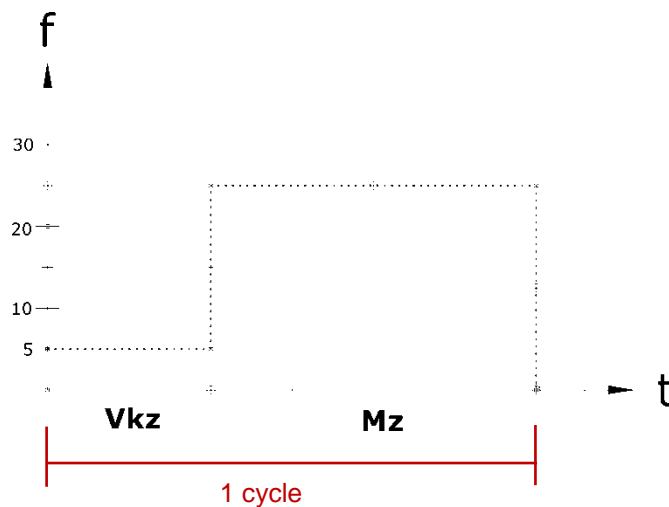


Fig. 30: Sequence of a single grinding cycle

The first grinding cycle consists of the pre-cooling and the set grinding time. All other grinding cycles consist of the intermediate cooling time and the set grinding time.

**Number of cryo cycles: n**

$$t_n = (Vkz + Mz) + [(n - 1)(Zkz + Mz)]$$

*t<sub>n</sub> = total grinding time*

$V_{zk}$  = precooling time

$M_z$  = grinding time

$Z_{kz}$  = intermediate cooling time

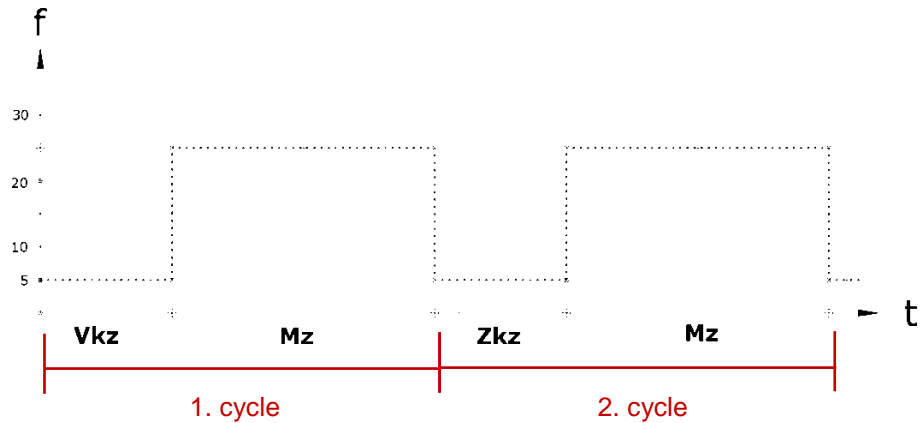


Fig. 31: Sequence of two grinding cycles

## 5.7 Programme Mode

### 5.7.1 Grinding without program


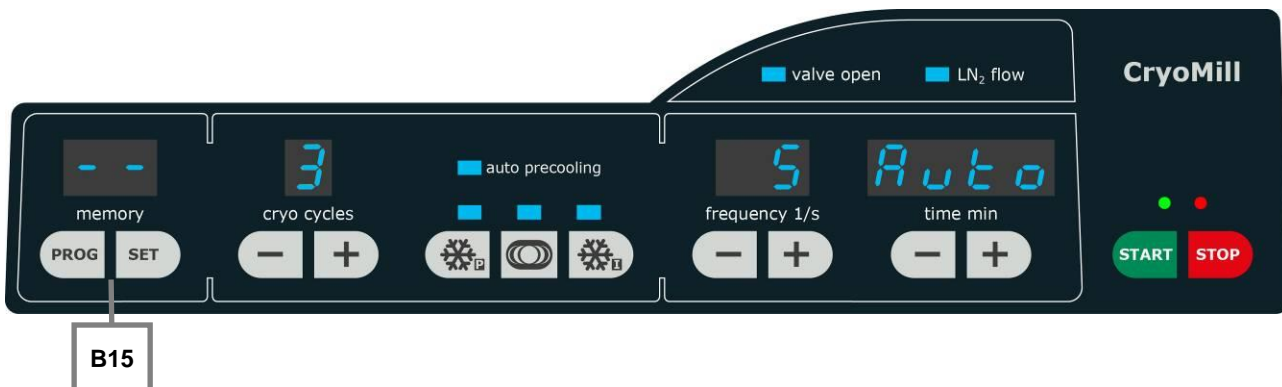


Fig. 32: Grinding without program

- Press the PROG button (**B15**) until “- -” is shown in the **memory** display.

The display advances one program each time the button is pressed

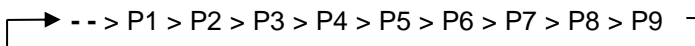

- The SET button (**B15**) is blocked in the “- -” mode

## 5.7.2 Grinding with program

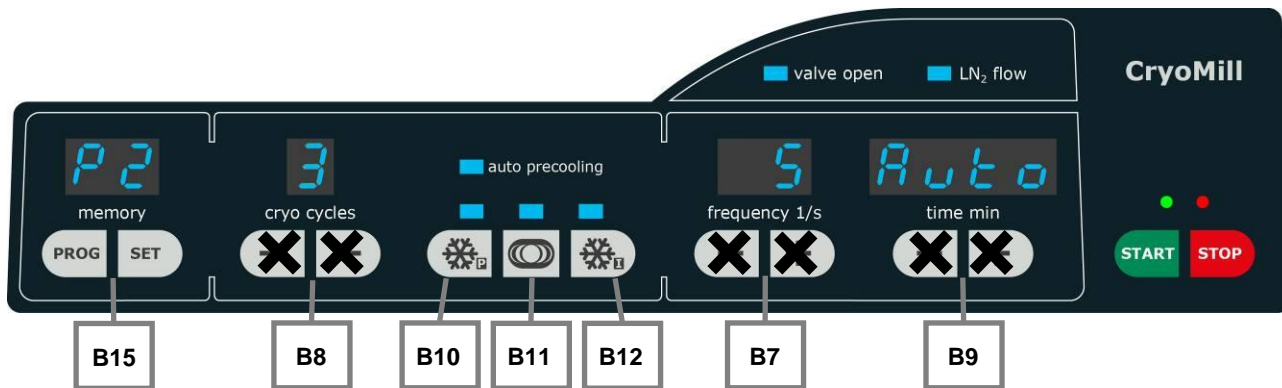


Fig. 33: Grinding with program

### 5.7.2.1 Selection of a program

- Press the PROG button (**B15**) until the required program is shown in the **memory** display.
- The display advances one program each time the button is pressed .  

→ - - > P1 > P2 > P3 > P4 > P5 > P6 > P7 > P8 > P9
- The buttons **B8**, **B7**, **B9** are blocked due to the selection of a program.

### 5.7.2.2 Display of the program grinding parameters

#### Precooling

- Press the precooling button (**B10**).
- The set time for precooling appears in the **time min (B9)** display.

#### Grinding time

- Press the grinding time button (**B11**).
- The set time for the grinding time appears in the **time min (B9)** display.

#### Intermediate cooling

- Press the intermediate cooling button (**B12**).
- The set time for the intermediate cooling appears in the **time min (B9)** display.

### 5.7.2.3 Create/change a program

- Press the PROG button (**B15**) until the required program (memory space) is displayed in the **memory** display.
- Press the SET button (**B15**).
- All displays and the corresponding LEDs for which settings are possible light up or flash.
- Set the parameters for grinding as described in the following chapters:
  - Number of grinding cycles
  - Precooling
  - Grinding time
  - Intermediate cooling time
  - Grinding frequency
- It is possible to cancel at any time using the PROG button (**B15**)



- Press the SET button (**B15**).
- Pressing the SET BUTTON (**B15**) again stores all parameters in the program.

#### 5.7.2.4 Delete program

- Press the PROG button (**B15**) until the program to be deleted is displayed in the **memory** display.
- Press the shortcut SET (**B15**) and - (**B8**) for at least 5 seconds.

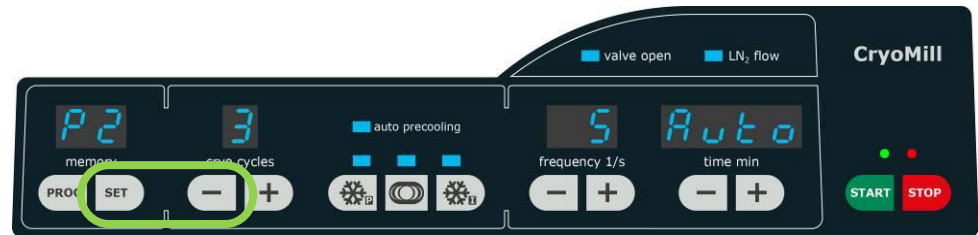


Fig. 34: Deleting a program

- The program is deleted and the parameters are assigned with the values “0” or “-”.
- It is not possible to start a deleted program.

## 5.8 Setting grinding cycles



Fig. 35: Setting the grinding cycles

- Switch the CryoMill on at the main switch.
- After switching on, the last used grinding parameters are displayed and can be used.
- You can set 0 (“-”) to 9 grinding cycles

#### NOTICE

If 0 (“-”) CryoCycles are set, the grinding is conducted without nitrogen. In this case the buttons **B10**, **B11** and **B12** are blocked.

- The first cycle consists of the precooling phase and grinding time
- Cycles 2 to 9 each comprise an intermediate cooling and grinding.
- Press the respective buttons (B8) to set the grinding cycles.
- **+** Pressing briefly increases the number.
- **-** Pressing briefly reduces the number.
- The display runs more quickly when pressed longer.

## 5.9 Setting the precooling time

### NOTICE

The machine runs at a fixed frequency of 5 Hz during the precooling and intermediate cooling. A range of from 30 seconds to 99 minutes can be set for precooling.

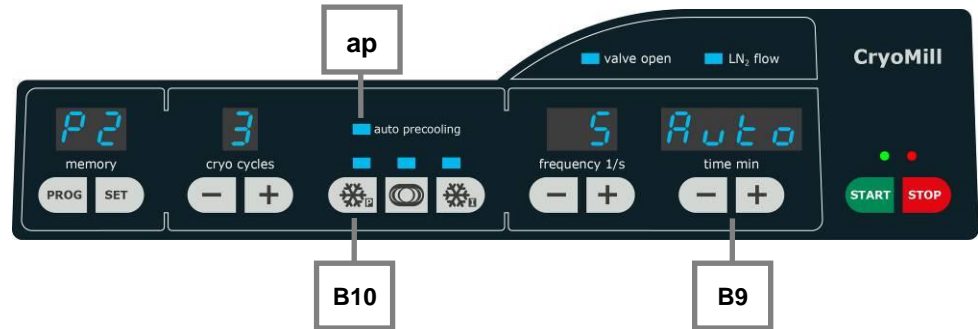


Fig. 36: Setting precooling

After switching on, the last used grinding parameters are displayed and can be used.

- Switch the CryoMill on at the main switch.
- Press the button (B10) to set precooling
  - The precooling LED lights up.

You can choose between two precooling versions:

### 5.9.1 Precooling with definable precooling time

- Use the buttons (B9) to set precooling.
- **+** Pressing briefly increases the duration.
- **-** Pressing briefly shortens the duration.

The display runs more quickly when pressed longer.

### 5.9.2 Automatic precooling

In the case of automatic precooling the grinding only starts once the entire system has been sufficiently cooled. This is checked by a sensor (see chapter >>Connecting the cooling agent supply <<)

#### Switch automatic precooling on

- Press the button (B10) for at least 2 seconds.
  - The auto precooling (ap) LED is lit up.
  - “Auto” is shown in the **time min** display.

#### Switch automatic precooling off

- Press the button (B10) for at least 2 seconds.
  - The auto precooling (ap) LED is extinguished.
  - “08:00” is shown in the **time min** display.

### 5.10 Setting the grinding time

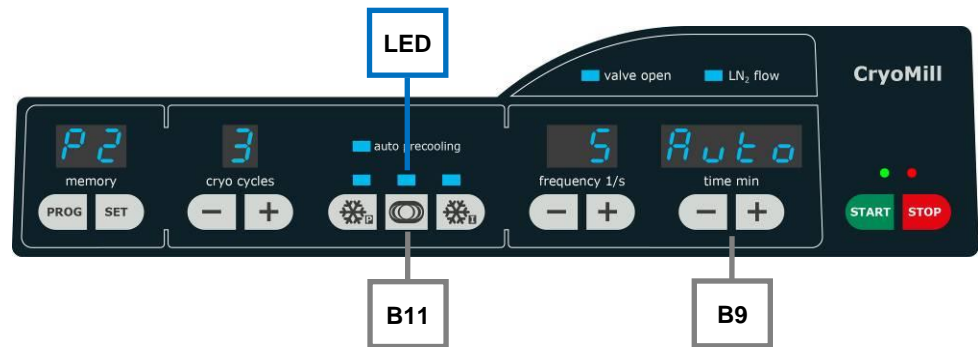


Fig. 37: Setting the grinding time

- Press the button (B11) to set the grinding time. The grinding time LED lights up.
- Use the buttons (B9) to set the grinding time
  - + Pressing briefly increases the duration.
  - Pressing briefly shortens the duration.
 The display runs more quickly when pressed longer.

*NOTICE*

In the case of grinding without cooling, the grinding time can be set directly using the buttons (B9).

### 5.11 Setting the duration of intermediate cooling

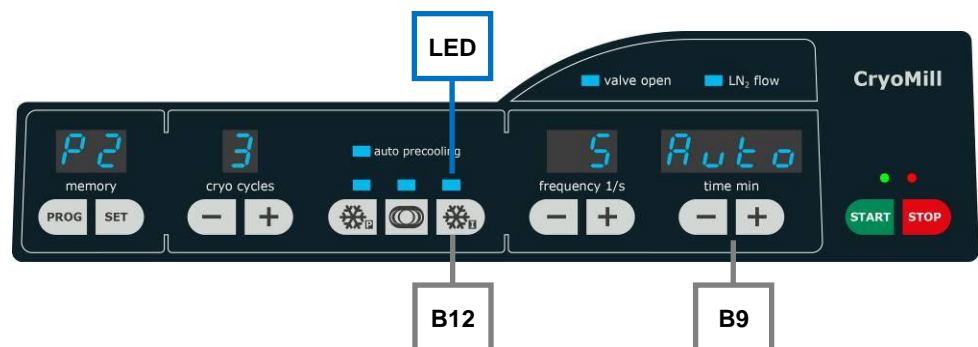


Fig. 38: Setting the intermediate cooling

- Press the button (B12) to set the intermediate cooling . The intermediate cooling time LED lights up.
- Use the buttons (B9) to set the intermediate cooling time (30 seconds to 99:00 minutes).
  - + Pressing briefly increases the duration.

- ▬ Pressing briefly shortens the duration.
- The display runs more quickly when pressed longer.

### 5.12 Setting the grinding frequency

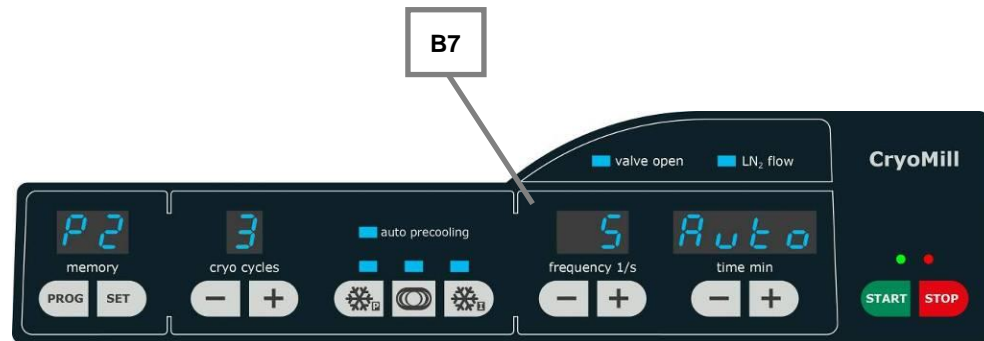


Fig. 39: Frequency 1/s

- Switch the CryoMill on at the main switch.
- It is only possible to adjust the grinding frequency if the button (B11) for setting the grinding parameters has been pressed and the LED (B11) is lit up.
- Press the button (B11).
- Press the respective buttons (B7) to set the grinding frequency .
- The selectable frequency range is from 5 to 30 Hz in steps of 1 Hz.
- + Pressing briefly increases the frequency.
- - Pressing briefly reduces the frequency.

The display runs more quickly when pressed longer.

The frequency of precooling and intermediate cooling is fixed at 5 Hz. The frequency cannot be adjusted while setting the precooling or intermediate cooling time.

### 5.13 Starting the grinding process



**Danger of injury and frostbite**

Escape of liquid nitrogen

- Liquid nitrogen has a temperature of -196 °C and may cause injuries similar to burns on skin or eye contact or frostbite.
- **Never switch the machine on without tightly closed cooling casing.**

*NOTICE*

The grinding process can only be started if the hood is closed.



Fig. 40: Starting the grinding process

- Press the START button to start the grinding process.
- The Start LED (green) lights up. At the same time the solenoid valve opens and liquid nitrogen can flow into the cooling system.

## 5.14 Interrupting and continuing the grinding process

### 5.14.1 Pausing grinding

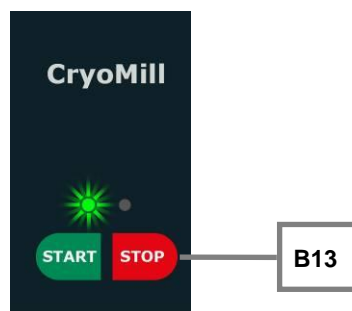


Fig. 41: Pausing grinding

- Press the STOP button (**B13**) once.
- Grinding is interrupted and the LED above the START button flashes.

### 5.14.2 Cancelling grinding

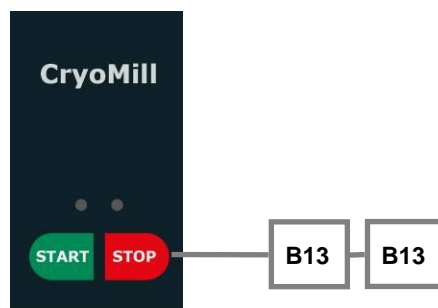


Fig. 42

- Press the STOP button (**B13**) twice.
- The grinding is ended and the two LEDs above the START and STOP buttons are off.

## 5.15 Stopping the grinding process

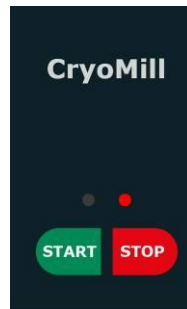


Fig. 43: Ending the grinding process

- Press the STOP button to end the grinding process.

The STOP LED (red) lights up. At the same time the solenoid valve closes.

In the event that the supply of liquid nitrogen has been interrupted during operation, this is indicated by the flashing of the remaining time LED (**B6**) at the end of the grinding time.

Pressing once interrupts the grinding process, for example to appraise the grinding sample. The remaining grinding time is still visible on the display.

By pressing the start button again, the mill continues running until the grinding time has completely expired.

Pressing the stop button twice aborts the grinding process. The machine is now in start mode.

By pressing the start button, the display is activated again and the grinding time is reset to the last start setting.

The duration can be changed during the grinding process. The grinding process is automatically ended on expiry of the grinding time. The display is reset to the last value started.

*NOTICE* If the milling process is aborted due to too little liquid nitrogen in the cooling circuit, the error E85 is displayed and the LED on top of the STOP button flashes. The error can be acknowledged by pressing the STOP button. The LEDs on top of the START and STOP buttons are alternately flashing and the remaining process time is displayed.

- To continue the grinding process (to do so enough liquid nitrogen must be present again), press the START button.
- To abort the grinding process and to display the original parameters again, press the STOP button once again.

An exception is the interruption during the automatic precooling stage (first cooling stage), before the process time (3 min) is counted down. As to this point in time, no process time is displayed, the time display (remaining process time) remains off in case of error E85.

## 5.16 Operating hours display

- Switch the device off at the mains switch.

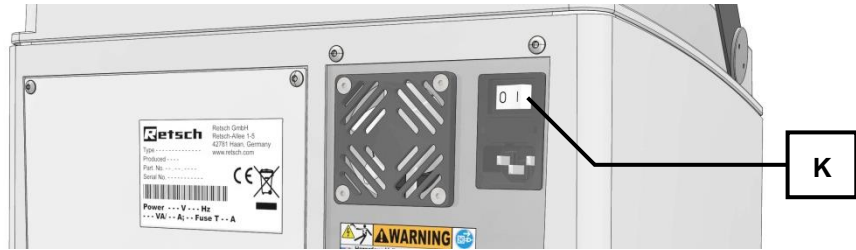


Fig. 44: Switching the device on and off

- Keep the buttons cryo cycles (**B8**) + and frequency (**B7**) - pressed simultaneously and with the buttons pressed, switch the device on at the mains switch.

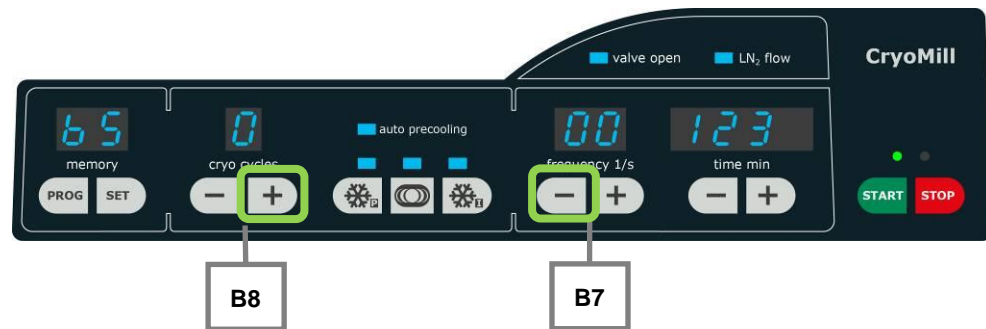


Fig. 45: Display of the operating hours

The total operating time in hours is displayed. The maximum value that can be depicted is 999999 hours. The letters bS are displayed in the first two places.

Examples:

bS0 00 012 = 12 hours operating time

bS1 23 456 = 123456 hours operating time

- Press the STOP button to exit the operating hours display

### 5.17 Operating software display

- Switch the device off at the mains switch.

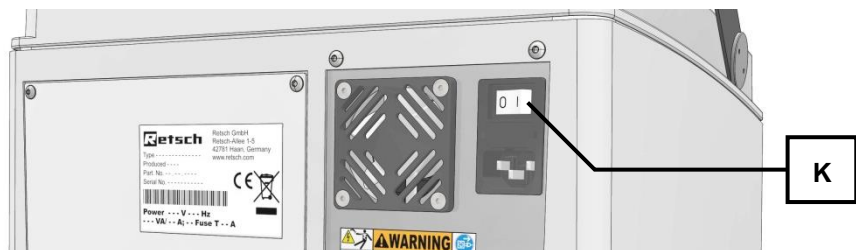


Fig. 46: Switching the device on and off

- Keep the buttons cryo cycles (**B8**) + and frequency (**B7**) + pressed simultaneously and...
- ...with the buttons pressed, switch on the device at the mains switch.

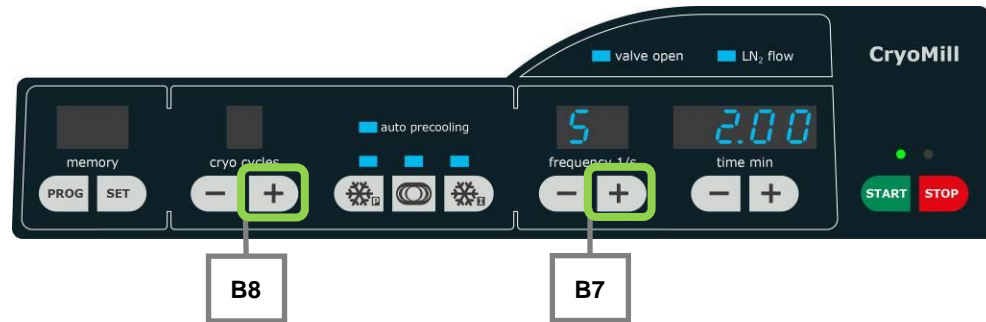


Fig. 47: Display of the software version

The current version of the operating software is displayed across the two right-hand display windows. The letter S is shown in first place.

Examples:

S 1.23 = Version 1.23 of the operating software

S 2.00 = Version 2.00 of the operating software

- Press the STOP button to exit the display of the operating software.

### 5.18 Replacing the machine fuses

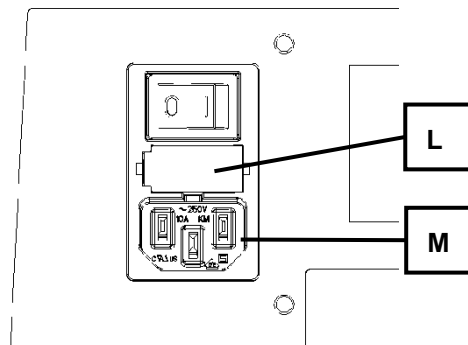


Fig. 48: Replacing the device fuses

The following glass fuses are needed for the CryoMill:

2 T4A

- Pull the mains plug from the connector (**M**).
- Press in the catch at the side. This releases the fuse holder (**L**) so that it can be pulled out.
- Replace the fuses.
- Push the fuse holder (**L**) in until it engages.



## 6 Working instructions

### 6.1 General

The CryoMill is an extremely modern, efficient product from Retsch GmbH.

Due to the large selection of accessories, the CryoMill has a wide range of potential applications in industrial and research laboratories.

It is primarily deployed in the chemical and pharmaceutical area, in mineralogy and biology etc.

---

 **CAUTION**

2.V0004

**Risk of explosion or fire**

Changing sample properties

- Consider that the properties and therefore also the hazardousness of your sample can change during the grinding process.
- **Do not use any substances in this device which carry the risk of explosion or fire.**

---

 **CAUTION****Danger of personal injury**

Dangerous nature of the sample

- Depending on the dangerous nature of your sample, take the necessary measures to rule out any danger to persons.
- **Observe the safety guidelines and datasheets of your sample material.**



### 6.2 The grinding process

- Insert your grinding jar filled with grinding sample and grinding ball/s into the grinding jar support as described in the chapter "Inserting the grinding jar".
- Close the hood.
- Set the grinding parameters.
- Open your external liquid nitrogen supply line.
- Start the device.

Do not open the hood during the grinding process.

The grinding jars are brought to a halt immediately by the built-in brake, however the grinding process can no longer be continued with the remaining grinding time.

The device must be restarted. The original parameters are available here.

### 6.3 Quantities of sample material and feed sizes

In addition to the instrument settings, the filling level of the grinding jar is also of crucial importance for a successful grinding process in ball mills of Retsch GmbH. **When grinding bulk materials, the grinding jar filling should consist of approximately one third of sample and one third of ball quantity.** The remaining third is the free grinding jar volume, which is required for the movement of the balls.

If an increase or decrease in sample volume is to be expected during the grinding process, the amount of sample can be adjusted within the range listed in the table. Thus, e.g. for voluminous materials such as wool, leaves, grasses and similar, a material filling level of 70 – 80 % is necessary. For wet grinding with grinding balls < 3 mm, the ball charge should be 60 % of the grinding jar volume.

Grinding jar nominal volume [ml]	Sample quantity [ml]	Max. feed size [mm]	Recommended ball filling (unit)						
			Ø5 mm	Ø7 mm	Ø10 mm	Ø12 mm	Ø15 mm	Ø20 mm	Ø25 mm
1.5	0.2 – 0.5	1	1 – 2						
5.0	0.5 – 2.0	2		1 – 2	-	-	-	-	-
10.0	2.0 – 4.0	4		5 – 7	1 – 2	1 – 2		-	-
25.0	4.0 – 10.0	6			5 – 6	2 – 4	1 – 2	-	-
35.0	6.0 – 15.0	6			6 – 9	4 – 6	2 – 3	1	
50.0	8.0 – 20.0	8			12 – 14	6 – 8	3 – 4	1	1

Safe-Lock reaction vessel	Sample quantity	Max. feed size	Recommended ball filling (unit)	
			Ø3 mm	Ø5 mm
1.5 ml	0.5 ml	1 mm	10 - 12	2 - 3
2 ml	0.7 ml	1.2 mm	14 - 16	3 - 4

### 6.4 Working instructions

The grinding jar and its content are cooled by means of liquid nitrogen (LN<sub>2</sub>) until it is possible to grind soft materials such as rubber.

Precooling is required until the ideal grinding temperature is reached. The drive operates at a lower frequency (5Hz) here to prevent the moving parts of the device from freezing. The fan is switched on during precooling and blows the air heated by the motor in the front through the openings of the rocker arms into the grinding chamber. This warm air prevents the rocker arm

bearing from cooling down. It is therefore also necessary to allow the fan to run on for a while after switching the grinding off.

In the preliminary tests an average precooling of approx. 10 to 15 minutes was determined. The grinding time was approx. 4 – 5 minutes.

The valve is switched off when the optimal temperature is reached. Grinding can begin as from this time and the valve will be switched on and off according to the temperature.

Grinding can also take place without cooling, in which case the cryo cycles are set to -.



Manual cooling is started by pressing the shortcut (B10) and (START B13). The valve remains open as long as the buttons are kept pressed.

## 7 Safety functions and fault display

### 7.1 Error Messages

Error code	(FEHLER) BESCHREIBUNG	DEFECT DESCRIPTION TRANSLATION
E10	ANTRIEB ÜBERLASTET	DRIVE OVERLOAD
E21	FEHLER DREHZAHL	FAILURE ROTATION SPEED
E22	FEHLER TASTATUR	FAILURE KEYPAD
E24	FEHLER VENTIL	FAILURE VALVE
E41	FEHLER DREHZAHLSENSOR	FAILURE SPEED SENSOR
E42	FEHLER TEMPERATURSENSOR 1	FAILURE TEMPERATURE SENSOR 1
E43	FEHLER TEMPERATURSENSOR 2	FAILURE TEMPERATURE SENSOR 2
E51	FEHLER SICHERHEITSSCHALTER	SATETY SWITCH DEFECTIVE
E85	KEIN STICKSTOFF	NO LIQUID NITROGEN
E86	LECKAGE	LEAKAGE
H40	MASCHINE SCHLIESSEN	CLOSE MACHINE

## 8 Cleaning, Wear and Maintenance

### 8.1 Cleaning

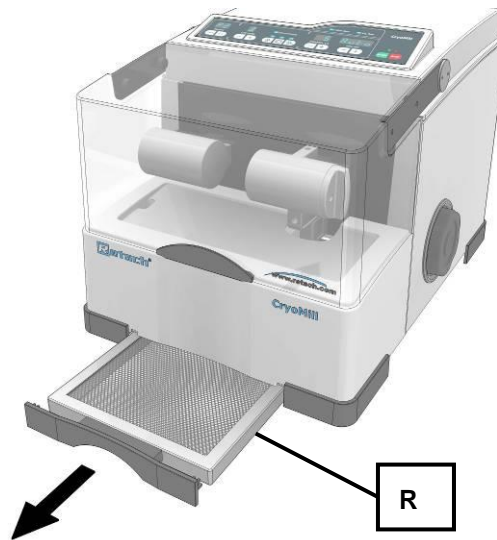


Fig. 49: Pulling out the condensation collecting filter

The condensation collecting filter must be cleaned regularly depending on the ambient conditions, the degree of use and the general air humidity.

- Pull out the condensation collecting filter (R).
- Clean the condensation collecting filter under running water with a little standard washing up liquid.
- Allow the collecting filter to dry and then insert in the device again.

### 8.2 Wear

The grinding tools may become worn, depending on the frequency of the grinding operation and the sample material. The grinding jars and the balls should be regularly checked for wear and replaced if necessary.

### 8.3 Wearing parts

---

 **CAUTION**

**Personal injury**  
Incorrect repairs

- This operating manual does not include instructions for repair.
  - **For your own safety repairs should be carried out only by Retsch GmbH or an authorised representative (service technicians).**
- 

### 8.4 Maintenance

The grinding tools may become worn, depending on the frequency of the grinding operation and the sample material. The grinding jars and the balls should be regularly checked for wear and replaced if necessary.

## 8.5 Checks

Correct functioning of the hood switch must be checked regularly.

- Switch on the machine with switch (K).
- Start grinding operation with START button.
- Lift the hood.

If there is an opening gap of a few cm the machine switches itself off and H40 appears in the display.

- Cancel H40 with the STOP button.



If this switching-off function does not take place the CryoMill must be checked immediately by Retsch-Service.

## 9 Safety and occupational health instructions

### CAUTION

#### Danger of injury to eyes and skin

Frostbites through liquid nitrogen

- Liquid nitrogen has a temperature of  $-196\text{ °C}$  and may cause injuries similar to burns on skin or eye contact or cause frostbite.
- **Always use goggles and wear protective gloves when opening the cooling casing and the grinding jar.**



### CAUTION

#### Use of liquid nitrogen

- Retsch GmbH rules out any liability claims that may arise in the use of liquid nitrogen or similar cooling agents.
- **Please observe the safety regulations of the cooling liquid supplier.**



### 9.1 General

This section summarises the general occupational health and safety guidelines for the handling of liquid nitrogen.

### CAUTION

#### Handling liquid nitrogen

General hazardous situations

- The following hazardous situations may arise in the handling of liquid nitrogen: oxygen deficiency situations, cryogenic burns, danger of explosion, oxygen enrichment
- **All users must be informed about the dangers of working with liquid nitrogen in order to be able to work safely.**

### 9.2 Safety instructions on the handling of liquid nitrogen

#### 9.3 Oxygen deficiency

The main components of air according to volumes are as follows:

- Oxygen O<sub>2</sub> 21 %
- Nitrogen N<sub>2</sub> 78 %
- Argon Ar 1 %

The gases contained in the atmosphere are not toxic. However, a change in concentration (in particular changes to the oxygen concentration) may have effects on life and burns. It is therefore essential for the air breathed in to contain sufficient oxygen (> 19 %).

Humans cannot detect changes in the composition of air within the time that will actually be necessary because the components are colourless and odourless.

### **9.3.1 Dangers**

A danger of suffocation exists as a result of the normal evaporation of the liquid nitrogen that pushes out the oxygen in the air. An oxygen deficiency is dangerous and can cause death through suffocation. The reaction of the organism to oxygen deficiency will differ greatly depending on the individual. It is not possible to provide precise and generally applicable information on the symptoms of an oxygen deficiency.

Example: under normal conditions (20° C; 1013 mbars) 1 l liquid nitrogen evaporates to produce 680 l nitrogen gas.

### **9.3.2 Causes**

An oxygen deficiency may arise during the following work or in the following conditions:

- Nitrogen as liquid or gas
- Natural evaporation of liquid nitrogen
- Refilling of liquid nitrogen
- Leaks in containers for liquid or gaseous nitrogen
- Defect in the air feed or outlet
- Tipping over of the container

This list is not complete.

### **9.3.3 Recommendations**

In order to prevent the danger of an oxygen deficiency, the following measures must be taken.

The vessel:

- must be kept in a vertical position.
- must be provided with a suitable insulating lid.
- must be protected from direct sunlight and may not be set up near heat sources.
- may not be transported in vehicles in filled state.
- must be protected from impact, knocks and rapid movements.
  
- Ventilate all installation rooms constantly and appropriately.
- Wear protective clothing (suitable gloves, goggles or face protection and safety shoes).
- Check the oxygen content of the room constantly.
- Always carry an oxygen metre.
- Only trained personnel may work with liquid nitrogen.

This list is not complete.



### 9.3.4 General conduct in the case of an accident

In the case of an accident from oxygen deficiency the following regulations should be observed.

- Secure the surrounding area to avoid any subsequent accidents.
- Act quickly.
- The rescuers must take measures to protect themselves (respiratory protection device).
- Move the injured persons from the hazardous area.
- Observe the company instructions for emergencies.
- Ventilate the rooms affected sufficiently.
- Investigate the cause of the accident.

This list is not complete.

## 9.4 Cryogenic burns

Liquid nitrogen is very cold (-196° C).

The surfaces of vessels that were in contact in liquid nitrogen (in particular during the filling process) may cause skin burns on contact.

### 9.4.1 Dangers

Cryogenic liquids may:

- bring about burns to the human body
- make specific materials (metal and plastic) that are not suitable for low temperatures brittle
- generate strong misting depending on atmospheric humidity

### 9.4.2 Causes

There are two types of cryogenic burns:

#### 9.4.2.1 Burns through splashes

When handling samples and in general when handling liquid nitrogen, personnel must protect themselves from splashes. They can cause cryogenic burns with serious consequential damage, in particular to eyes and face.

#### 9.4.2.2 Burns through contact

Contact of the skin with cold material causes frostbite or cryogenic burns.

The interior of vessels or the samples may never be touched or held with bare hands.

### 9.4.3 Recommendations

In order to prevent the danger of burning the following points must be observed:

- Never bring cryogenic liquids in contact with the skin.
- Never touch the cold non-isolated or iced walls of a vessel.
- Wear personal protective equipment (suitable gloves, goggles or face protection and safety shoes).
- Keep the vessel upright.

- Use suitable material (e.g. metal hose or PTFE hose) to refill.
- Train personnel.

This list is not complete.

#### **9.4.4 General rules of conduct for splashes with liquid nitrogen**

##### **9.4.4.1 To the eyes**

- Rinse eyes with much water for 15 min.
- Follow the company instructions for emergencies.
- Consult a doctor

##### **9.4.4.2 On the skin**

- Do not rub.
- If possible remove or loosen clothing.
- Slowly and gradually warm the parts affected.
- Do not apply anything to the burnt area.
- Follow the company instructions for emergencies.
- Consult a doctor.

Both lists are not complete.

### **9.5 Danger of explosion**

#### **9.5.1 Dangers**

The evaporation of liquid nitrogen may lead to an overpressure in the vessel.

#### **9.5.2 Causes**

An increase in the vessel may be attributable to:

- incorrect set-up (use of a closable lid)
- icing on the neck and on the insulating lid

This list is not complete.

#### **9.5.3 Recommendations**

To avoid the danger of explosion:

- Always use suitable insulating lids (pay attention to waste gas opening).
- Observe filling levels to avoid the formation of ice on the insulating lid.
- Set up the vessel in dry and roofed rooms.
- Monitor atmospheric humidity in the installation room.
- Check vessel regularly for the collection of condensed water.
- Check vessel regularly for damage to surface or material damage

This list is not complete.

#### 9.5.4 General conduct in the case of an accident

In the case of an accident from oxygen deficiency the following regulations should be observed.

- Secure the surrounding area to avoid any subsequent accidents.
- Act quickly.
- The rescuers must take measures to protect themselves (respiratory protection device).
- Move the injured persons from the hazardous area.
- Observe the company instructions for emergencies.
- Ventilate the rooms affected sufficiently.
- Investigate the cause of the accident.

This list is not complete.

### 9.6 Oxygen enrichment

#### 9.6.1 Dangers

Oxygen enrichment may increase the danger of explosion and fire.

#### 9.6.2 Causes

Oxygen may be condensed from the air if liquid nitrogen is used and similarly be liquefied because the boiling point of oxygen (approx.  $-183^{\circ}\text{C}$ ) is above that of nitrogen ( $-196^{\circ}\text{C}$ ).

#### 9.6.3 Recommendations

The following points must be considered in the case of oxygen enrichment:

- Do not smoke.
- Where possible, keep easily inflammable materials away from the vessel.
- Remove all sources of fire (naked flame and light, spark producers, matches, lighters etc.).
- Ventilate installation rooms constantly and appropriately.
- Clean floor regularly.
- Train personnel.
- Wear personal protective equipment.
- Check oxygen content constantly.
- Always carry an oxygen metre.

This list is not complete

## **9.7 Surrounding area of the machine**

### **9.7.1 Rooms**

The room in which the machine is located must:

- permit operation without danger to the employees
- have a constantly running and appropriate ventilation system
- have a level and non-porous floor that is also able to bear the load of the vessel
- have safety datasheets on liquid nitrogen visible for all to see
- prevent unauthorised access
- permit the safe filling of the vessel
- permit accessibility to the vessel for inspection, cleaning and maintenance

This list is not complete.

## 10 Disposal

Please observe the respective statutory requirements with respect to disposal.

Information on disposal of electrical and electronic machines in the European Community.

Within the European Community the disposal of electrically operated devices is regulated by national provisions that are based on the EU Directive 2002/96/EC on Waste Electrical and Electronic Equipment (WEEE).

Accordingly, all machines supplied after 13.08.2005 in the business-to-business area to which this product is classified, may no longer be disposed of with municipal or household waste. To document this they have the following label:

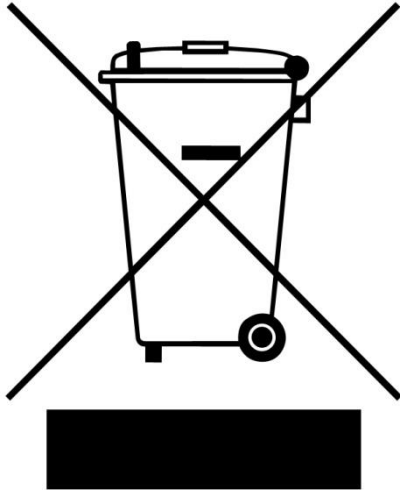


Fig. 50: Disposal label

Since the disposal regulations within the EU may differ from country to country we would request you to consult your supplier.

## 11 Index

<b>4</b>	
4 A.....	48
<b>5</b>	
5ml grinding jar position.....	36
<b>A</b>	
Ambient temperature.....	14
Amperage.....	20
Application-technical information.....	7
Assign/change memory space.....	40
Atmospheric humidity.....	15
Automatic precooling.....	42
<b>B</b>	
Bar code.....	20
bS.....	47
<b>C</b>	
Capacity.....	20
Causes.....	56, 57, 58, 59
CE marking.....	20
Changes.....	7
Checks.....	54
Clean.....	53
Cleaning.....	53
Cold air outlet.....	25
Condensation.....	24
Condensation collecting filter.....	17, 53
Confirmation.....	11
Connecting the coolant feed.....	21
Connection cable.....	20
Cooling jacket.....	33
Cooling jacket seal.....	22
Copyright.....	7
Create/change a program.....	40
cryo cycles.....	38
Cryogenic burns.....	57
<b>D</b>	
Danger of explosion.....	58
Dangers.....	56, 57, 58, 59
Degree of protection.....	13
Delete program.....	41
Depth.....	13
Device designation.....	20
Diagram of the cooling system.....	23
Dimensions.....	13
Dismantling the grinding jar retaining device.....	32
Display of the grinding parameters.....	40
Display of the operating hours.....	47
Display of the software version.....	48
Displays.....	30
Disposal.....	61
Disposal label.....	20
Disposal label.....	61
<b>E</b>	
Electrical connection.....	20
Emission value with regard to workplace.....	12
Emissions.....	12
Ending the grinding process.....	46
Eppendorf.....	36
Eppendorf reaction vial.....	37
Error messages.....	52
Exhaust hose.....	26
Explanations of the grinding cycles.....	37
Explanations of the safety instructions.....	8
External fuse.....	20
<b>F</b>	
Fill level.....	27
Final fineness levels.....	27
Flashing of the remaining time LED.....	46
frequency 1/s.....	38
Frequency range.....	44
Front view.....	28
Function of the hood switch.....	54
Functions.....	30
Fuse holder.....	48
Fuse strength.....	20
Fuse type.....	20
<b>G</b>	
G 1/4 inch.....	24
General.....	49, 55
General conduct in the case of an accident.....	57, 59
General rules of conduct for splashes with liquid nitrogen.....	58
General safety instructions.....	9
Glass fuses.....	48
Grinding cycles.....	41
Grinding frequency.....	44
Grinding jar	
inserting.....	31
Grinding jar retaining device.....	32
Grinding jar volumes.....	27

Grinding time .....	43	<b>O</b>	
Grinding time .....	38	Opening gap .....	54
Grinding with cooling .....	38	Operating controls.....	30
Grinding with program .....	40	Operating elements and displays .....	28
Grinding without cooling .....	37	Operating hours .....	47
Grinding without program .....	39	Operating hours display .....	46
Grip .....	19	Operating hours display .....	47
<b>H</b>		<b>Operating software display</b> .....	47
H40 .....	54	Operating the device.....	27
Height.....	13	Operating time .....	47
Hood .....	49	Operation .....	16
<b>I</b>		Optimal grinding jar filling .....	27
Insert the adapter for 4/2 x 5 ml grinding jar.....	35	Original parameters .....	49
Inserting a 5ml grinding jar .....	35, 36	Outlet for gaseous nitrogen .....	25
<b>Inserting and removing CryoMill special grinding jar</b> .....	32	Oxygen deficiency.....	55
Inserting reaction vials .....	37	Oxygen enrichment.....	59
Inserting the cooling jacket seal .....	35	<b>P</b>	
Inserting the cooling jacket sealing plug.....	35	Packaging .....	14
Inserting the grinding jar .....	33	Part number .....	20
Inserting the grinding jar adapter .....	36	Power version .....	20
Installation.....	14	Precooling .....	42
Installation height.....	15	Precooling .....	38
Installation of the device .....	15	Precooling with definable precooling time .....	42
Installation site		pressure relief .....	25
conditions .....	14	pressure relieve .....	25
Intermediate cooling .....	43	Programme mode .....	39
Intermediate cooling time.....	38	property damage.....	8
Interrupting and continuing the grinding process	45	Protective equipment.....	12
<b>L</b>		Pulling out the condensation collecting filter.....	17
Laboratory mill .....	12	<b>Q</b>	
Layer of ice .....	24	<b>Quantities of sample material and feed sizes</b>	50
Liquid nitrogen supply interrupted .....	46	<b>R</b>	
$L_{pAeq}$ .....	12	Rated power.....	13
<b>M</b>		reaction vial.....	36
Mains frequency .....	20	Rear view .....	28
Maintenance .....	53	Recommendations .....	56, 57, 58, 59
Manual .....	7	Recovery of the samples .....	27
Manufacturer's address .....	20	Regulations for the place of installation.....	20
Materials .....	27	Remaining grinding time .....	49
Maximum relative humidity .....	15	Remove the grinding jar.....	34
Moderate or mild injury .....	8	Removing the liquid nitrogen feed line.....	24
<b>N</b>		Removing the liquid nitrogen feed line.....	24
Name of machine model.....	27	Repair .....	10
Noise values .....	12	Repair instructions .....	7
Notes on the Manual.....	7	Replacing the machine fuses.....	48
Number of fuses .....	20	Required floor space.....	13
		Revision status.....	7
		Rooms.....	60

**S**

Safety and occupational health instructions .....	55
Safety functions and fault display .....	52
Safety instructions on the handling of liquid nitrogen .....	55
Safety warnings .....	8
Seal .....	22
Securing the grinding jar .....	33
Selection of a program .....	40
Sequence of a single grinding cycle .....	38
Sequence of two grinding cycles .....	39
Serial number .....	20
serious injury .....	8
Service Address .....	10
<b>Setting grinding cycles</b> .....	41
Setting precooling .....	42
<b>Setting the duration of intermediate cooling</b> .....	43
Setting the grinding cycles .....	41
<b>Setting the grinding frequency</b> .....	44
Setting the grinding time .....	43
<b>Setting the precooling time</b> .....	42
Single grinding cycle .....	38
Software .....	48
Special grinding jars .....	27
Starting the grinding process .....	44
Starting the grinding process .....	45
Stopping the grinding process .....	46
Summary table of device parts .....	29
Supply line pressure .....	24
Surrounding of the machine .....	60

**T**

T 4 A .....	48
Target group .....	27
Target group .....	9
Technical data .....	12
Temperature fluctuation and condensed water ..	14
<b>The grinding process</b> .....	49
threaded connector .....	24
Total operating time .....	47
Transport .....	14, 16
Transportation lock .....	16
removing .....	16
Type plate .....	20
description .....	20
Type plate lettering .....	20

**U**

Use of the device for the intended purpose .....	27
Using the opening aid .....	34

**V**

vapour plume .....	15
Version of the operating software .....	48
Voluminous materials .....	27

**W**

Wear .....	53
Wearing parts .....	53
Weight .....	13
Width .....	13
Working instructions .....	49
Working instructions .....	50

**Y**

Year of production .....	20
--------------------------	----



## CRYOGENIC MIXER MILL

CryoMill | 20.749.xxxx

### EU DECLARATION OF CONFORMITY

Herewith we declare, represented by the signatory, that the above mentioned device complies with the following directives and harmonized standards:

#### **Machinery Directive 2006/42/EC**

Applied standards, in particular:

DIN EN ISO 12100      Safety of machinery

#### **EMC Directive 2014/30/EU**

Applied standards, in particular:

DIN EN 55011      Industrial, scientific and medical equipment - Radio-frequency disturbance characteristics - Limits and methods of measurement

DIN EN 61000-3-2      Electromagnetic compatibility (EMC)

DIN EN 61000-3-3      Electromagnetic compatibility (EMC)

DIN EN 61326-1      Electrical equipment for measurement, control and laboratory use - EMC requirements

#### **Low Voltage Directive 2014/35/EU**

Applied standards, in particular:

DIN EN 61010-1      Safety requirements for electrical equipment for measurement, control and laboratory use

#### **Authorized person for the compilation of technical documents:**

Dr. Loredana Di Labio (technical documentation)

Furthermore, we declare that the relevant technical documentation for the above mentioned device has been compiled according to Annex VII Part B of the Machinery Directive, and we undertake to submit this documentation on request to the market surveillance authorities.

**In case of a modification of the device not previously agreed with Retsch GmbH, as well as the use of unauthorised spare parts or accessories, this declaration will lose its validity.**

Retsch GmbH



Dr. Ing. Frank Janetta, Team Leader R&D Department

Haan, 05/2016









**Retsch**<sup>®</sup>

**Copyright**

© Copyright by  
Retsch GmbH  
Retsch-Allee 1-5  
42781 Haan  
Germany