

Instructions for Use

Charger shell 2 – pipette holder 2

English (EN)

1 About this manual

1.1 Notes on this manual


This document supplements the operating manual for the product. This document does not replace the operating manual.

1. Read these instructions for use before using the product.
2. Also read the operating manual for the device with which you are using the accessory.


The dates in this manual correspond to the international date format as specified in the ISO 8601 standard. All dates are shown in the format YYYY-MM-DD or YYYY-MM.

2 Installation

2.1 Attaching the charger shell or pipette holder to the wall

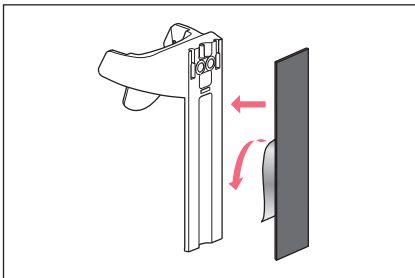

1. Clean the smooth glass or stainless steel surface (e.g., with Ethanol) and allow it to dry.
2. Remove the protective foil from one side of the adhesive tape.
3. Press the adhesive tape firmly onto the rear side of the charger shell or pipette holder.

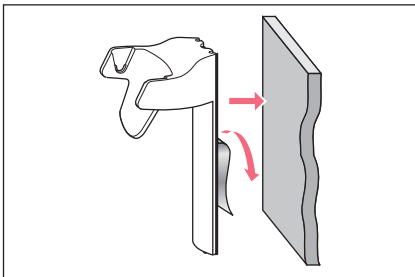

4. Remove the protective foil from the second side of the adhesive tape.
5. Press the charger shell or pipette holder firmly against the cleaned surface. Apply no load to the charger shell or pipette holder for at least 24 hours.

2.2 Removing the charger shell or pipette holder from the wall

1. Twist the charger shell or pipette holder and pull off the adhesive strip.

2 Instructions for Use

Charger shell 2 – pipette holder 2

English (EN)

2.3 Inserting the charger shell or pipette holder in the carousel

1. The charger shell or pipette holder is inserted into the carousel without adhesive strip.

3 Operation

3.1 Charging the electronic dosing device

The charge function is not carried out via the pipette holder with wall mounting. The dosing device is charged via its corresponding power supply unit.

4 Maintenance

4.1 Cleaning the charger shell or pipette holder

Material:

- Cloth
- Soap-based cleaning agent

1. Clean any contamination with a damp cloth and mild cleaning agent.

5 Technical data

5.1 Material

Open surfaces	Acrylonitrile styrene acrylate (ASA/PC)
Charging contacts	Nickel-plated coating

6 Ordering information

Charger shell

Description	Order no.
Charger Shell 2 for one Eppendorf Xplorer®, Xplorer plus or Xplorer plus Move It® electronic pipette for Charger Carousel 2 or Charger Stand 2, with charging function	3116 602 007
for one Multipette® E3/E3x or Multipette® stream/Xstream elec- tronic multi-dispenser for Charger Carousel 2 or Charger Stand 2, with charging function	3116 603 003

Pipette holder

Description	Order no.
Pipette Holder 2 white for Eppendorf Research® 3 pipettes for Pipette Carousel 2, Charger Carousel 2, Stand 2 or wall mounting, sticky tape included, without charging function	3116 000 295

Description	Order no.
Pipette Holder 2 for one Eppendorf Research®, Eppendorf Research® plus, Eppendorf Reference®, Eppendorf Reference® 2, Biomaster® pipette or epMo- tion® Tip Tool for Pipette Carousel 2, Charger Carousel 2, Stand 2 or wall mounting, sticky tape included, without charging function	3116 000 112
for one Eppendorf Xplorer®, Xplorer plus or Xplorer plus Move It® electronic pipette for Pipette Carousel 2, Charger Carousel 2, Stand 2 or wall mounting, sticky tape included, without charging function	3116 000 120
for one Multipette® E3/E3x or Multipette® stream/Xstream elec- tronic multi-dispenser for Pipette Carousel 2, Charger Carousel 2, Stand 2 or wall mounting, sticky tape included, without charging function	3116 000 139
for one Multipette® M4 multi-dispenser for Pipette Carousel 2, Charger Carousel 2, Stand 2 or wall mounting, sticky tape included, without charging function	3116 000 147