

Gage Block Assemblies

Parr gage block assemblies combine the function of an:

- Inlet / Outlet Valve
- Pressure Gage
- Safety Rupture Disc

into a compact assembly which can be attached to the head of any pressure vessel with a single connecting tube. There is a threaded socket in the block for a gas connection with a pressure hose or tubing using a type "A" coned pressure fitting. Coned fittings are also used on the thick-walled tube which connects the block to the pressure vessel. The valve in this assembly controls the gas flow into the vessel and the gage shows the pressure in the vessel when the valve is closed. A rupture disc with a burst rating matched to the gage is installed in the block.

Two Sizes

These gage block assemblies are made in two styles which differ only in the size of the pressure gage and the size of the tube connection to the

pressure vessel. The smaller 4316 Gage Block Assembly is usually furnished with a 3-1/2 inch diameter pressure gage and includes a 3 inch long connecting tube with type "A" coned fitting. This assembly is normally installed on pressure vessels up to 600 mL.

The 4317 Gage Block Assembly has a 4-1/2 inch diameter gage and type "B" coned fitting on a 3 inch long connecting tube. This assembly is normally installed on vessels of 1000 mL and larger.

The gage block, connecting tube and pressure gage on all of these assemblies are normally made of type 316 stainless steel, but they can also be made of Alloy 400 on special order. Longer tubes can be furnished on special order.

Ordering Information

Part No.	Description
4316	with Type "A" Connector
4317	with Type "B" Connector


Gages for Parr Pressure Vessels

Pressure Range		4-1/2" Diameter	3-1/2" Diameter
psi	bar	Gage Number	Gage Number
0-100	0-7	56HCPA	593HCP1AD
0-200	0-14	56HCPB	593HCP2AD
0-600	0-40	56HCPC	593HCP6AD
0-1000	0-70	56HCPD	593HCPD
0-2000	0-140	56HCPF	593HCPF
0-3000	0-210	56HCPG	593HCPG
0-5000	0-350	56HCPH	593HCP50AD
0-7500	0-517	56HCP75AD	NA
0-10000	0-700	56HCPK	NA

When ordering, specify the catalog number, pressure gage diameter and range.

Available Gages

Any of the gages shown in this table can be installed on these gage block assemblies.

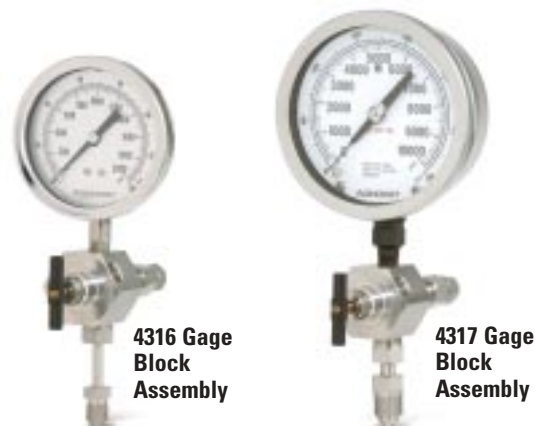

4316 Gage Block Assembly

4317 Gage Block Assembly