


INTRODUCING NEW

PRISMTM mini

Labnet International is proud to introduce our new personal mini centrifuge, the *PRISM MINI*. Reliable and easy-to-use, the *Prism Mini* is compact and economical, allowing each workstation to be equipped with a "personal" centrifuge. The *Prism Mini* can be used for a wide range of molecular biology separations and quick spins. The centrifuge was designed to be extremely compact with an exceptionally small footprint. It includes 2 quick release interchangeable rotors, which spin up to 6000 rpm's. The electronic brake and self opening lid allow for quick deceleration to reduce handling time. The *Prism Mini's* maintenance-free drive is remarkably quiet and the unique design allows for excellent air-flow, to protect temperature sensitive samples.

- Quick release rotor exchange
- Increased sample capacity
- Rapid sample processing
- One-step sample retrieval

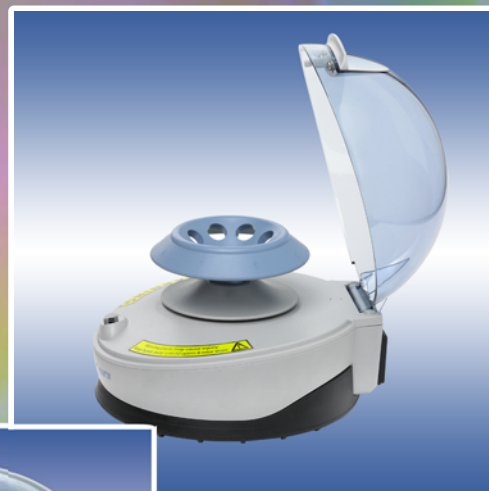

Specifications:

Capacity:	8 place tube rotor & 4 place PCR Strip rotor
Rotor Speed:	6000 rpm
RCF Force:	2000 x g
Dimensions:	5.2 x 5.2 x 4.4 inch 13.2 x 13.2 x 11.2 cm
Weight:	1.2lbs / .55 kg
Electrical:	120V or 230V, 50/60Hz

Ordering Information:

Catalog No.	Description
C1801	Prism Mini Centrifuge, 120V
C1801-230V*	Prism Mini Centrifuge, 230V
C1205-8	Adapters for 0.5 mL tubes, pk 8
C1206-8	Adapters for 0.4 mL tubes, pk 8
C1222-8	Adapters for 0.2 mL tubes, pk 8
*Specify cord type by adding -EU or -UK to the end of the catalog number	