

Package rotary-vane pump RZ 6 +FO +VS 16 +VACUU·VIEW extended

[Link to the website](#)

Product description

These powerful rotary vane pumps feature a particularly compact design and low weight for pumps of this capacity. They are the ideal solution for a wide range of laboratory and process applications that require low ultimate vacuum at medium to increased gas flow rate. The two-stage rotary vane pump RZ 6 is also available as a package including pump, oil mist filter FO, valve VS 16 and vacuum gauge VACUU·VIEW extended.

Performance features

- high flow rates even at vacuum levels approaching ultimate vacuum
- high water vapor tolerance due to efficient gas ballast; very good ultimate vacuum even with gas ballast
- vacuum-tight at switch-off; external anti-suckback valve not needed
- large oil volume: Long intervals between oil changes
- ease of maintenance due to telescopic design

Ordering data

Product name	Nominal mains voltage mains frequency	Plug type	Articleno.
RZ 6 +FO +VS 16 +VACUU-VIEW extended	230 V 50-60 Hz	CEE	20798160

Items supplied

Pump oil filled and completely mounted, ready for use, with manual, with oil mist filter FO, valve VS 16 and vacuum gauge VACUU-VIEW extended.

Accessories

- Rubber vacuum tubing DN 10 mm (20686002)
- Metal vacuum tubing with small flanges KF DN 16 (20673336)
- Separator inlet side AK for RE/RZ 5/6/2.5 (20698006)
- Maintenance kit RE 2.5 - RZ 6 (20648977)


Technical data

Max. pumping speed 50 Hz (m³/h)	5.7 m³/h
Max. pumping speed 60 Hz (cfm)	4.0 cfm
Ultimate vacuum	2×10^{-3} mbar / 1.5×10^{-3} torr
Ultimate vacuum with gas ballast	1×10^{-2} mbar / 0.75×10^{-2} torr
Ultimate partial vacuum without gas ballast	4×10^{-4} mbar / 3×10^{-4} torr
Water vapor tolerance	40 mbar / 30 torr
Number of stages	2
Vacuum controller	VACUU·VIEW extended
Permissible ambient temperature range (operation)	12 - 40 °C
Permissible ambient temperature range (storage)	-10 - 60 °C
Max. outlet pressure (abs.)	1.1 bar
Oil capacity min./max.	0.34 / 0.73 l
Inlet connection	Small flange KF DN 16 / hose nozzle DN 10 mm
Outlet connection	Hose nozzle DN 8-10 mm
Rated motor power	0.18 kW
Rated speed 50/60 Hz	1500/1800 min ⁻¹
Protection class IEC 60529	IP 40
Dimensions (L x W x H)	420 mm x 230 mm x 240 mm
Weight	17.8 kg
Noise (sound pressure level) in dBA at 50 Hz or 1500 rpm/62% (VARIO)/1500 rpm (VARIO-SP)/12500 Upm (VACUU·PURE®)	50 dBA
NRTL certification	No

Pumping speed graph at 50 Hz



Pump down graph at 50 Hz (10 l volume)

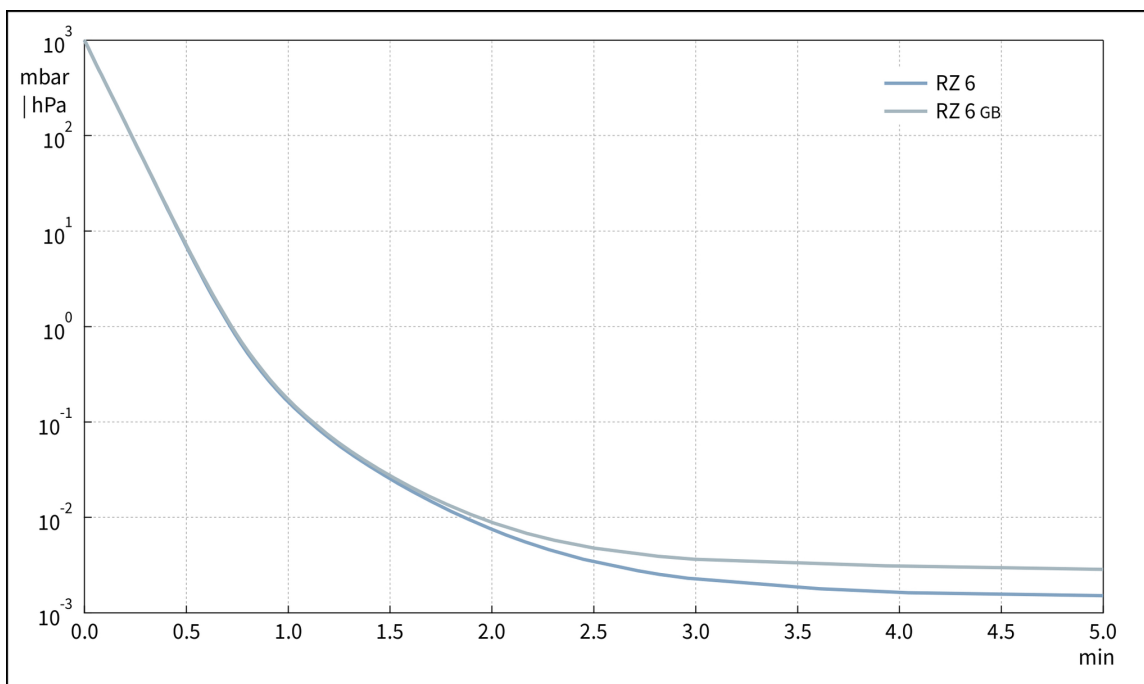

Pumping speed graph at 60 Hz



Pump down graph at 60 Hz (10 l volume)



Dimension sheet



VACUUBRAND GMBH + CO KG
P.O. Box 1664
97877 Wertheim
T +49 9342 808 5550
F +49 9342 808 5555

info@vacuubrand.com
www.vacuubrand.com

VACUUBRAND®

VACUUBRAND is part of the Brand group, a manufacturer of high-quality and innovative laboratory equipment and vacuum technology with a global team of approx. 1,000 employees.

With highly motivated and qualified employees, we manufacture the world's most comprehensive product family for vacuum generation, measurement and control for the rough and fine vacuum range in the laboratory.

All with one goal:

Enabling our customers to run reliable, safe and efficient processes. This is why we commit ourselves to continuous innovation, first-class performance, and outstanding quality.

Engineering excellence is central to our ability to deliver products and solutions renowned for their quality, durability, and reliability. Innovative technology, high quality materials and extensive testing deliver unmatched performance - engineered and manufactured in Germany.

Together with you, our experts will find the best solution for your application and provide you with first-class service along the way.

© VACUUBRAND GMBH + CO KG

BRAND (Shanghai) Trading Co., Ltd.
Shanghai, China

T +86 21 6422 2318
info@brand.com.cn
china.brand.com.cn

BRAND Scientific Equipment Pvt. Ltd.
Mumbai, India

T +91 22 42957790
customersupport@brand.co.in
www.brand.co.in

BRANDTECH® Scientific, Inc.
Essex, CT. United States of America

T +1 860 767 2562
info@brandtech.com
www.brandtech.com