


Data Sheet

Mixing / Shakers

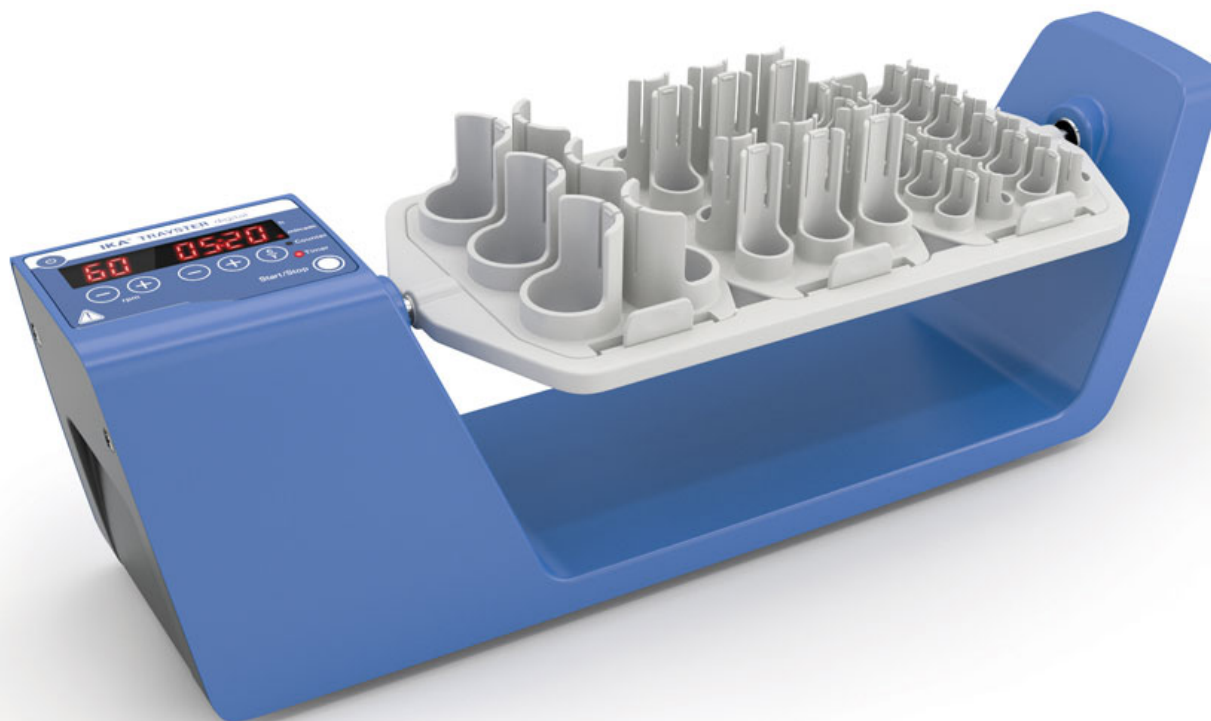

Trayster digital

Digital overhead shaker with vertical rotation for smooth but effective mixing with infinitely adjustable speed.

Ideal for biological samples, particularly blood samples. Suitable for separating and mixing processes of powdered and liquid samples in Eppendorf tubes and Greiner tubes up to 50 ml. The easy removable attachments allow simple and quick cleaning of the unit of spilled samples.

- Timer: Countdown, adjustable from 1 s to 99h 59min
- Counter: Display of shaking time
- Different attachments allow a variety of applications
- Usable at 4 to 50 °C in incubators
- Suitable for continuous operation
- Error code display
- Easy operation by foil keypad
- Speed range: 5 - 80 rpm

Accessories: TR 1 Tube insert, TR 2 Tube insert, TR 3 Tube insert, TR 4 Tube insert

Technical Data	
Type of movement	overhead
Permissible shaking weight (incl. attachment) [kg]	2
Permissible ON time [%]	100
Speed min (adjustable) [rpm]	5
Speed range [rpm]	0 - 80
Speed display	7 segment LED
Speed control	1 RPM steps
Timer	yes
Timer display	7 segment LED
Time setting min. [s]	1
Time setting max. [min]	5999
Operating mode	timer and continuous operation
Dimensions (W x H x D) [mm]	460 x 180 x 140
Weight [kg]	1.6
Permissible ambient temperature [°C]	4 - 50
Permissible relative humidity [%]	80
Protection class according to DIN EN 60529	IP 21
Voltage [V]	100 - 240
Frequency [Hz]	50/60
DC Voltage [V=]	24
Current consumption [mA]	1000
Ident. No.	0004006000