

BIOCLEAN FUSION

Sterile Polychloroprene Gloves
S-BFAP

About

BioClean Fusion Sterile Polychloroprene Cleanroom Gloves contain no natural latex proteins, are comfortable and thin to reduce wearer fatigue and increase tactility. The antistatic glove is double-donnable and features a beaded cuff for strength.

Specifications

COMPATIBILITY:

ISO Class 4 & EU GMP Grade A

LENGTH:

290mm (11.4")

MATERIAL:

Polychloroprene

PROTEIN LEVEL:

Polychloroprene contains no natural latex proteins

SURFACE:

Finger-textured

SHAPE:

Ambidextrous

COLOUR:

Green

STERILIZATION:

Gamma irradiation

TACTILE

CHEMICAL RESISTANT*

DOUBLE-DONNABLE

Features

- Increased tactility
- Powder-free & latex-free
- Antistatic
- Beaded cuff for strength
- Non-particulating EasyTear packaging


EN 421:2010 ISO 374-5:2016



EN ISO 374



CE 0598

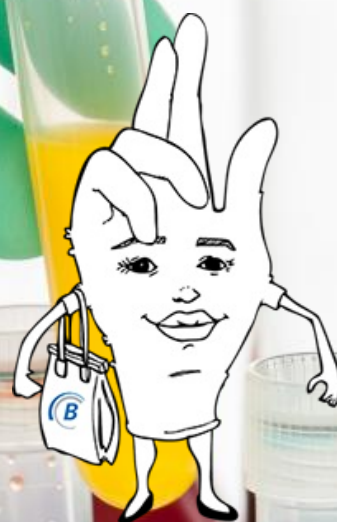

Ms.G.Love™ represents BioClean's range of cleanroom gloves



Quality Standards

- Manufactured in a facility holding ISO 9001:2015 and (EU)2016/425 Module D certifications
- Conforms to Category 3 - Complex Design Personal Protective Equipment PPE Regulation (EU)2016/425
- Processed in an NEBB certified ISO Class 4 environment
- Complies with the requirements of European Glove Standards EN420:2003 +A1:2009, EN374-1:2016, EN374-5:2016 and EN421:2010
- Physical properties comply with European Medical Glove Standard EN455-2:2015
- Sterilized in accordance with ISO 11137-1:2006
- Suitable for contact with all foodstuffs. Regulation (EC)1935/2004 and (EC) 2023/2006.
- Tested against permeation standard ASTM D6978-05 for cytotoxic drugs

Physical Properties

CHARACTERISTICS	VALUE	TEST METHOD
Freedom From Holes	AQL 0.65 Performance Level 3	EN 374-2:2015
Glove Length +/-10mm (mm/inches)	290/11.4	EN 420:2003 + A1:2009
Single Wall Thickness Minimum (mils/mm)	Cuff: 2.76/0.07 Palm: 3.94/0.10 Finger Tip: 4.72/0.12	EN 455-2:2015
Force At Break During Shelf Life	≥9 N	EN 455-2:2015

PALM WIDTHS (mm) ± 5mm						
SIZE	5.0-5.5	6.0-6.5	7.0-7.5	8.0-8.5	9.0	10.0
mm	76	86	95	106	115	120
TEST METHOD				EN 420:2003 + A1:2009		

Sterilization

Method: Gamma irradiation

Minimum Dose: 25kGy

Sterility Assurance Level: 10⁻⁶

Shelf Life & Storage

Five (5) years from date of manufacture.
Store in a dry, cool place (<40°C) away
from direct sunlight and fluorescent
light.

Glove Cleanliness Characteristics

PARTICLES					TYPICAL PARTICLE COUNT			TEST METHOD
≥ 0.5μm(counts/cm ²)					850			IENT-RP-CC005.4
TYPICAL EXTRACTABLE IONS								TEST METHOD
ANIONS	Fluoride	Chloride	Nitrite	Bromide	Nitrate	Phosphate	Sulphate	IENT-RP-CC005.4
(μg/cm ²)	ND	0.194	ND	ND	0.884	ND	0.025	
CATIONS	Lithium	Sodium	Ammonium	Potassium	Calcium	Magnesium	Zinc	
(μg/cm ²)	ND	0.099	0.003	0.315	0.652	0.007	ND	
ND = Not Detected, NT = Not Tested								

TO ORDER

Re-order code

S-BFAP-5055
S-BFAP-6065
S-BFAP-7075
S-BFAP-8085
S-BFAP-9090
S-BFAP-1010

Size

5.0-5.5
6.0-6.5
7.0-7.5
8.0-8.5
9.0
10.0

Packing

One pair per inner PE wallet; one wallet per sealed EasyTear PE pouch; 10 pouches per sealed outer PE bag; 20 outer bags per lined carton (200 pairs)

For additional information visit us at www.bioclean.com, or call us at +44 1638 663338

Nitritex Ltd

Minton Enterprise Park
Oaks Drive
Newmarket, Suffolk
CB8 7YY
United Kingdom

S-BFAP/1904/PDS28

Ansell

Neither this document nor any other statement made herein by or on behalf of Ansell should be construed as a warranty of merchantability or that any Ansell product is fit for a particular purpose. Ansell assumes no responsibility for the suitability or adequacy of an end user's selection of gloves for a specific application.

*Warning: No glove provides complete protection against all chemicals. Users must test the gloves against the particular chemicals and environment in which they will be used.

**Please see product validation pack or contact Ansell customer service for specific data on use of gloves with cytotoxic drugs. Gloves used for protection against such drugs must be selected specifically for the type of chemicals used.

Ansell, ® and TM are trademarks owned by Ansell Limited or one of its affiliates. © 2012 – 2019 Ansell Limited. All Rights Reserved.