


Application Number 126820
Method HPLC
Results

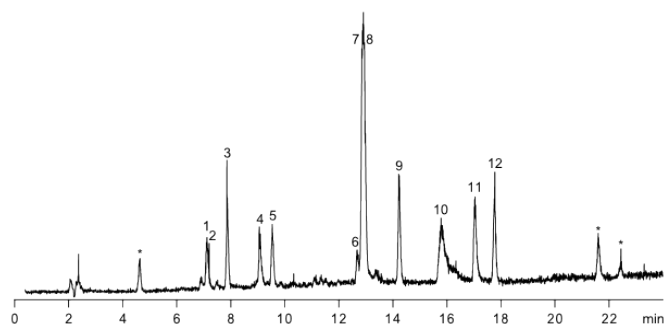

Legend:

Peaks (m/z):

- 1 - 403.8 (+2)
- 2 - 339.7 (+2)
- 3 - 317.7 (+2)
- 4 - 793.0 (+2)
- 5 - 728.9 (+2)
- 6 - 479.1 (+3)
- 7 - 390.2 (+3)
- 8 - 779.5 (+1)
- 9 - 482.9 (+2)
- 10 - 545.5 (+3)
- 11 - 713.5 (+3)
- 12 - 1006.0 (+3)
- * - unknown

Substances tryptic digest of cytochrome C (bos taurus)

Product(s)

Phase

REF

Webshop

NUCLEODUR
C18 ec

760186.20


[Shop now](#)

Matrix

Sample(s)	tryptic digest of cytochrome C
Conditions	Eluent A: 0.1 % formic acid in water Eluent B: 0.1 % formic acid in acetonitrile Gradient: 5 % B to 60 % B in 25 min Flow rate: 0.3 mL/min Temperature: 25 °C Injection volume: 5.0 µL
Detection	MS ESI positive Fullscan, 300 - 1010 m/z
Note	
Author	HPLC department
Source	MACHEREY-NAGEL, Germany 2014
Disclaimer	Concerning questions to literature applications please contact the cited author/s. Chemicals mentioned in applications are not manufactured or sold by MACHEREY-NAGEL. Applications mentioned in MACHEREY-NAGEL literature or databases are provided for informational purposes only. MACHEREY-NAGEL does not warrant that all applications have been tested in MACHEREY-NAGEL laboratories using MACHEREY-NAGEL products. MACHEREY-NAGEL does not warrant the correctness of the applications.

Share via Mail

 Print View

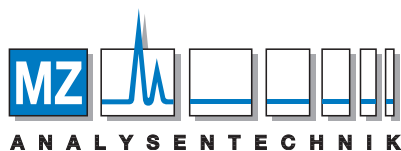

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de